

ELGA

MEDICA EDI - US
Operator Manual

Copyright Note

The information contained in this document is the property of VWS (UK) Ltd. and is supplied without liability for errors or omissions.

No part of this document may be reproduced or used except as authorized by contract or other written permission from VWS Ltd. The copyright and all restrictions on reproduction and use apply to all media in which this information may be placed.

VWS Ltd. pursue a policy of continual product improvement and reserve the right to alter without notice the specification, design, price or conditions of supply of any product or service.

© VWS (UK) Ltd. 2012

All rights reserved.

Publication ref: MANU40685

Version 2 – 10/15

ELGA LabWater is a trading name of VWS (UK) Ltd.

ELGA® and PURELAB® are registered trademarks.

MinnCare is a trademark of Minntech Corporation.

EfferSan™ is a trademark of Activon Inc

TABLE OF CONTENTS


1.	INTRODUCTION	4
	1.1 Use of this Manual.....	4
	1.2 Customer Support	4
	1.3 HEALTH AND SAFETY NOTES	4
	1.4 Electricity	4
	1.5 Pressure	4
	1.6 Hydrogen.....	5
	1.7 Sanitization Chemicals.....	5
	1.8 Control of Substances Hazardous to Health (COSHH).....	5
	1.9 Consumables.	5
2.	INSTALLATION AND COMMISSIONING	6
	2.1 Installation	6
	2.2 Commissioning.....	6
3.	CONTROLS	7
	3.1 User Recognition Keys.....	8
	3.2 System Preferences	8
	3.3 Setting Consumable Replacement Reminders	16
4.	OPERATION	18
	4.1 Night service/Operational day (screen 030).....	18
	4.2 Continuous Recirculation (24/7) (screen 029)	18
5.	MAINTENANCE	19
	5.1 Replacing Composite Vent Filter (L136)	19
	5.2 Replacing Purification Cartridge LC207 or Conditioning Cartridge LC206.....	20
	5.3 Replacing Ultra-Microfilter (LC109).....	21
	5.4 Replacing UV lamp (LC105)	22
	5.5 Replacing Protek L2 (LC177) or Carbon Filter (LC123) Pre-treatment.	23
	5.6 Replacing E- Cartridge (LC181).....	24
	5.7 RO modules	24
6.	SANITIZATION	25
	6.1 Liquid sanitization.....	25
	6.1.1 RO sanitization (Step 1) RO sanitization.....	25
	<i>Set Sanitization reminder</i>	27
	6.2 Tablet sanitization	28
	6.2.1 RO sanitization (Step 1).	28
	<i>Set Sanitization reminder</i>	30
7.	EMERGENCY BYPASS	31
	7.1 Emergency bypass.....	31
	7.2 Calculation for DI Cartridge (LC207) Cartridge life ...	32
8.	TROUBLE SHOOTING	33
9.	KEY TO CONTROL PANEL	34

10. TECHNICAL SPECIFICATION 35

11. WARRANTY/CONDITIONS OF SALE 37

12. USEFUL CONTACT DETAILS 39

1. INTRODUCTION

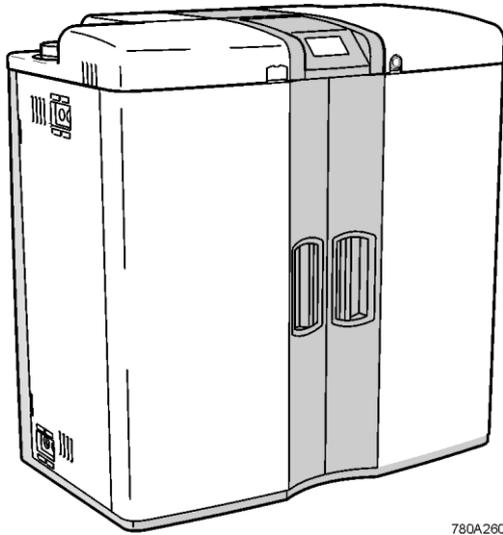

780A260

1.1 Use of this Manual

ELGA LabWater strives to produce manuals that are as simple and accurate as possible. However, should you feel they can be improved in any way please email us at info@elgalabwater.com

This manual guides you through the installation, commissioning and basic operation of the **MEDICA EDI - US** allowing you to obtain a guaranteed supply of purified water to meet your requirements. For more detailed information refer to the Operator Manual supplied on CD.

1.2 Customer Support

Service support and consumable items are available from your local supplier or distributor. Please contact info@elgalabwater.com for further information.

1.3 HEALTH AND SAFETY NOTES



WARNING! WARNINGS ARE GIVEN WHERE FAILING TO OBSERVE THE INSTRUCTIONS COULD RESULT IN INJURY OR FATALITY.



CAUTION! Cautions are given where failure to observe the instructions could result in damage to the equipment, associated equipment and processes.

1.4 Electricity

It is essential that the electrical supply to the **MEDICA EDI - US** is isolated before any items are changed or maintenance work performed.

The appliance coupler (mains lead) connected to the rear of the unit can be removed to isolate the power supply. If access to this is restricted then it is recommended that access to the supply socket is easily available to disconnect the electrical supply.

The unit must be earthed. Use of the recommended appliance coupler will ensure adequate earth protection.



WARNING! IF THE EQUIPMENT IS USED IN A MANNER NOT SPECIFIED BY ELGA, THE PROTECTION PROVIDED BY THE EQUIPMENT MAY BE IMPAIRED.

1.5 Pressure

The water supply should be isolated and residual pressure released prior to removal of any components or carrying out work on the unit.

Switch off the process and relieve pressure in the distribution loop by opening a point of use or dispense tap.

1.6 Hydrogen

Small quantities of Hydrogen gas are present in the drain line. Drain must be in a ventilated area.

1.7 Sanitization Chemicals

During the automatic sanitization cycles, EfferSan™ multi-purpose disinfecting tablets and Minncare cold sterilant are used and relevant safety guidance is included in this Manual. Please refer to the manufacturer for material safety data sheets.

EfferSan™ and Minncare cold sterilant are EPA registered as a sterilant, high level disinfectant, and sanitizer.

EfferSan™ and Spent Minncare cold sterilant acidic and requires normal neutralization as specified by your local state and local regulations.

1% of Minncare cold sterilant has a pH of 3.5.

1.8 Control of Substances Hazardous to Health (COSHH)



Material safety data sheets covering the various replaceable components are available upon request. Contact your local ELGA LabWater distributor.

1.9 Consumables.

Consumables	Typical Service Life*	Max. Shelf Life
LC105 (UV lamp)	12 months	5 years
LC109 (UMF)	6 months	2 years
LC136 M2 (Vent filter)	6 months	2 years
LC207 (DI Cartridge)	3-12 months	2 years
LC206 (Conditioning Cartridge)	2-6 months	2 years
LC 177 (Protek L2)*	12 months	2 years
LC193 (RO cartridge 120 L/hr)	3 years	2 years
LC181 (E-Cartridge)	3 years	2 years

* Service Life is an estimate only, and will depend on the application and feed water quality. Care should be taken to order the correct consumable items.

* Restricted use based on feed water conditions. Refer to technical specification on page 31.

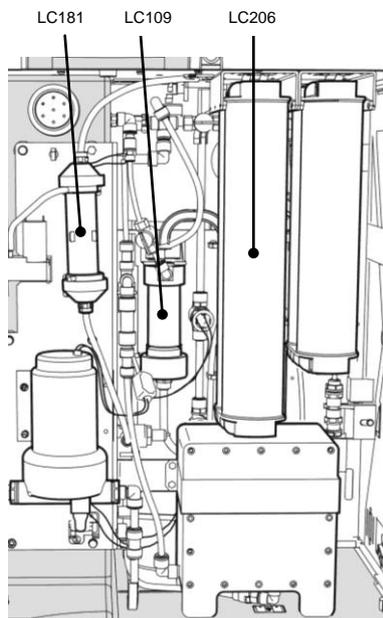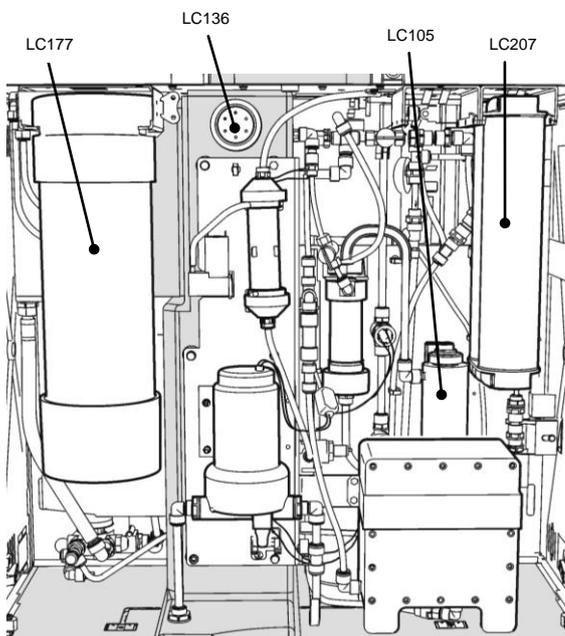

2. INSTALLATION AND COMMISSIONING

2.1 Installation


The unit is supplied with a quick start manual that allows the unit to be quickly installed and commissioned.

It is recommended that a trained ELGA representative carry out the installation of the product.

2.2 Commissioning

The unit is supplied with the software set in a commissioning mode that must be completed before being used to supply the application.

It is recommended that a trained ELGA representative complete the installation of the product.

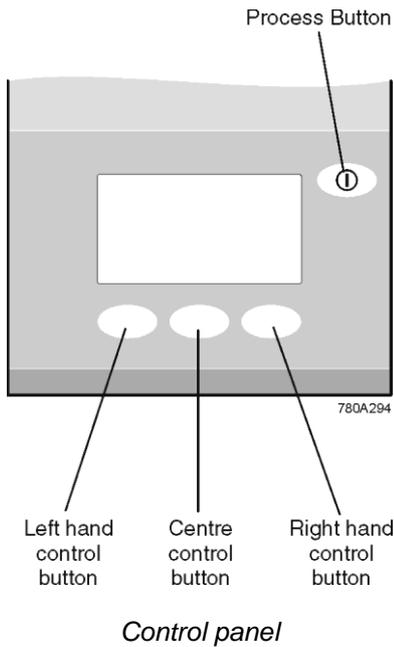

3. CONTROLS

The **MEDICA EDI - US** control panel has a range of control icons. General icons are as follows. Further icons are described in the appropriate sections a complete listing is included in Section 9.

BUTTON	ICON	Function
PROCESS		Turns the unit on/off
LEFT		Menu
		Scroll
		Shift
		Up
CENTRE		Accept
		Replacement dates
RIGHT		Mute Alarm
		Down

3.1 User Recognition Keys


Your unit is supplied with the following User keys:

Master Key	(Black)	1off
User Key	(Blue)	4off
Sanitization Key	(Green)	2off

The Master Key should be stored in a safe place. The Master Key controls the access level of the other User Keys.

User Keys only have access to customer preference screens.

Sanitization Keys start the sanitization process and stop the general user from accidentally starting sanitization.

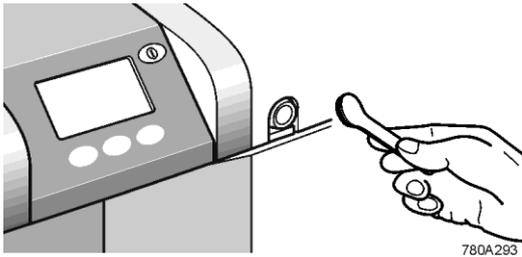

Present the Master or User Key



CAUTION! If the Master Key is lost a new Key can only be programmed by ELGA Service Engineers.

3.2 System Preferences

When the **MEDICA EDI - US** unit is started for the first time after completing the commissioning routine the following steps should be carried out to set up your choices:

Note: Your choices can be changed and implemented during normal operation. It is not necessary to stop the unit.

Step 1 - System access (User Key)

The User Key prevents unauthorized access to specific settings. This ensures consistent system performance and operation.

Note: The User Key does not prevent access to the PROCESS function in case of emergency.

1. SWITCH on the main electrical supply to start the controller set-up sequence. This takes several seconds.

PRESENT the 'Master Key' (black) or the 'User key' (blue) to the reader ensuring clean contact of both metallic components.


The display will show a padlock 

REMOVE 'Key' from the reader.

The display will show a key 

PRESS  button.


Note: If no buttons are pressed the system will relock after 5 seconds

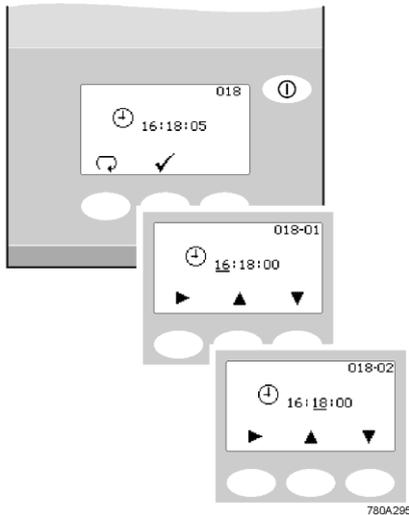

Setting clock time

Step 2 - Clock (screen 018)

1. PRESS  to edit time
OR
PRESS  to proceed to step 3.
2. PRESS and HOLD  to increase or  to decrease hour.
3. PRESS  to step cursor onto minute.
4. PRESS  to increase or  to decrease minutes.
5. PRESS .

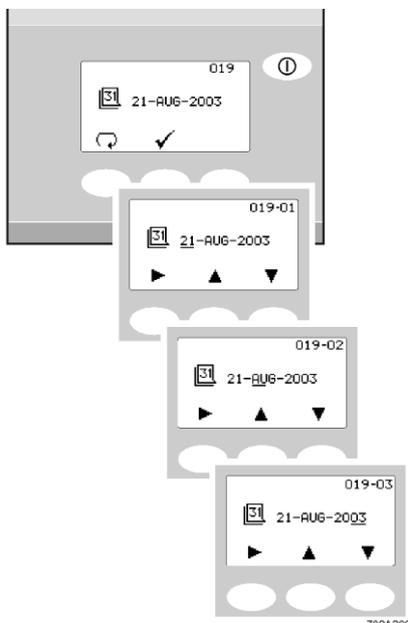


Set date screens

Step 3 - Date (screen 019)

The date is used to initiate change reminders, the date will appear on printed records.

1. PRESS  to edit date
OR
PRESS  to proceed to step 4.
2. PRESS and HOLD  to increase or  to decrease day.
3. PRESS  to step cursor onto month.
4. PRESS  to increase or  to decrease month.
5. PRESS  to step cursor onto year.
6. PRESS  to increase or  to decrease year.
7. PRESS .

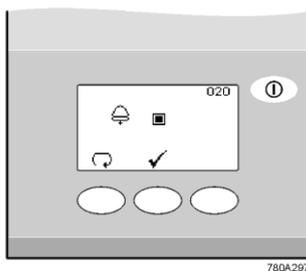


Audible alarm enable screen

Step 4 - Audible alarm enabled/disabled (screen 020)

This screen provides the option of either turning on the audible alarm, causing it to sound (while the alarm icon flashes) or turning off the audible alarm meaning it will remain quiet

1. PRESS  to change mode ( = ON)
OR
PRESS  to proceed to step 5.
2. PRESS .

Note: The visual alarm cannot be disabled.

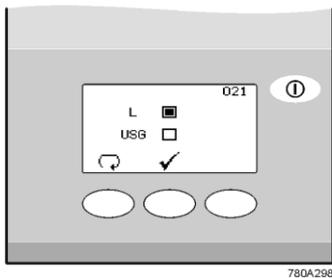


Set units of measure - volume

Step 5 –Set units of measure - Volume (screen 021)

This screen allows units of water volume to be set, to either **Litres** or **US Gallons**. This only indicates the volume in the reservoir.

1. PRESS to change (L or USG)
OR
PRESS to proceed to step 6.
2. PRESS .

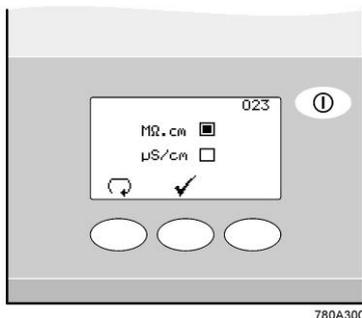


Water purity unit setting

Step 7 - Water purity unit setting (screen 023)

This screen allows preferred displayed units of water purity to be set, to either, **MΩ.cm** or **µS/cm**. This only affects the quality measurement in the recirculation loop.

1. PRESS SCROLL to change mode (MΩ.cm or µS/cm)
OR
PRESS to proceed to step 8.
2. PRESS .

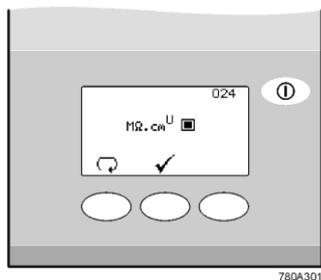

Uncompensated water quality screen

Step 8 - Uncompensated water quality (screen 024)


A U will indicate uncompensated readings (recirculation loop only) in the normal process screen.

1. PRESS to change (■ = Uncompensated water quality ON)
OR
PRESS to proceed to step 9.
2. PRESS .

Note: This mode is for validation purposes only and should not be selected as standard.

Step 9 - RO water quality alarm settings QS1 (screen 025)

This screen is used to set the value at which the RO water quality alarm will sound. This alarm does not stop RO production but can extend or start flush routines.

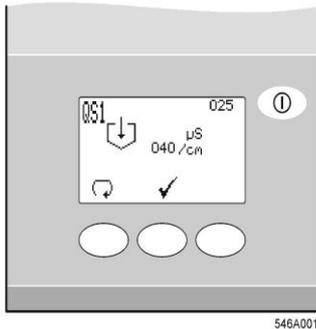


Alarm settings QS1 screen

1. PRESS  to select alarm point (Increments of 10, ranging from 20 to 100μS/cm). If '0' is selected the alarm will be turned 'OFF'.
OR
PRESS  to proceed to step 10.
2. PRESS .

Step 10 - Product water purity alarm settings QS2 (screen 026)

This screen is used to select the value at which the product water purity alarm will activate. The alarm does not stop the unit and will automatically reset if the purity level recovers

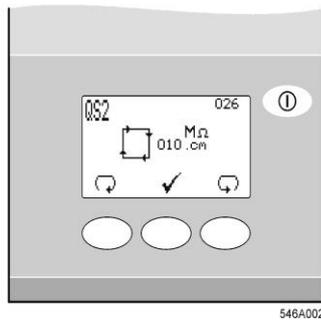


Alarm settings QS2 screen

1. PRESS  to select alarm point (increments of 1, ranging from 1 to 10 MΩ.cm).
OR
PRESS  to proceed to step 11.
2. PRESS .

Step 11 - RO water temperature alarm settings TS1 (screen 027)


This screen is used to select the value at which the RO water temperature alarm will activate. The alarm does not stop the unit and will automatically reset if the temperature returns below the set point.

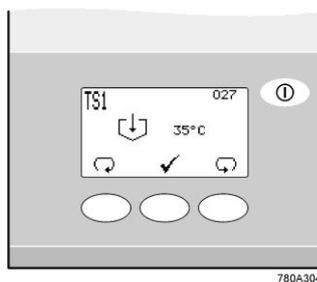

Alarm settings TS1 screen

1. PRESS  to select alarm point (increments of 1, ranging from 20 to 50°C)
OR
PRESS  to proceed to step 12.
2. PRESS .

Step 12 - Product water temperature alarm settings TS2 (screen 028)

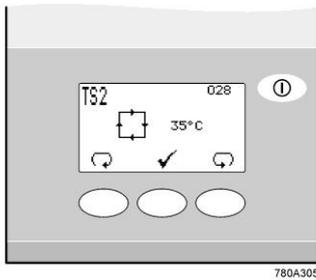


Alarm settings TS2 screen

This screen is used to select the value at which the product water temperature alarm will activate. The alarm does not stop the unit and will automatically reset if the temperature returns below the set point.

1. PRESS to select alarm point (increments of 1, ranging from 20 to 50°C)
OR
PRESS to proceed to step 13.
2. PRESS .

Note: If the water temperature rises above 55°C in the system it will alarm, stop and await operator intervention - alarm will be reset once power is removed for 5 seconds and then reinstated

Step 13 - Continuous operation (screen 029)

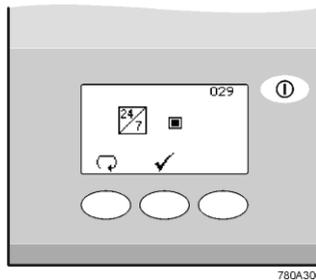

Continuous operation screen

Continuous operation may be required in particular circumstances or when demand for water fluctuates.




CAUTION! During long periods of continuous operation and low usage the water temperature will rise. It is recommended that this function is only used when water usage is on average >50 litres/hr and water is used every day.

1. PRESS to change (■ = ON)
OR
PRESS to proceed to step 14.

Note: Proceed to step 16 if continuous operation is selected (■ = ON).

Step 14 - Periods of operation (screen 030)

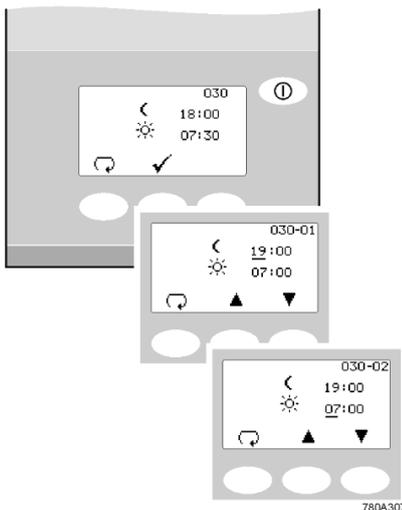

Setting period of operation

For efficiency and to reduce the likelihood of heat build up the normal operational hours can be selected.

During 'OFF' periods the unit will automatically re-circulate water for a period of 10 minutes every two hours this will maintain water purity within the distribution loop.


1. PRESS to edit night service start
OR
PRESS to proceed to step 15.
2. PRESS to increase or to decrease time in increments of 30 minutes.
3. PRESS to step to night service end.
4. PRESS to increase or to decrease time in increments of 30 minutes.
5. PRESS to enter times.
6. PRESS .

Note: Night service will only be operational if Continuous Operation is not enabled.

Step 15 - Operational days selection (screen 031)

Select the days that the **MEDICA EDI - US** is required to operate by highlighting the relevant box.

Monday = 1, Sunday = 7

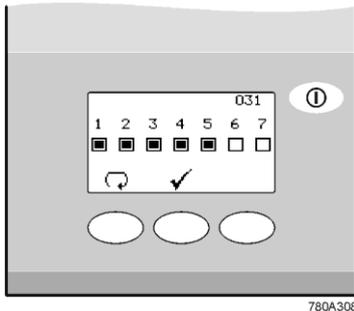


Day selection screen

1. PRESS to enter selection screen
OR
PRESS to proceed to step 14.
 2. PRESS to highlight box 1 (■ = Monday enabled)
OR
PRESS to step to box (2).
- Note: A highlighted box indicates that the unit will be operational on that day between the times set in step 14.*
3. REPEAT, step 15 - item 2, to select further operating days or PRESS until the appears.
 4. PRESS .

Note: During the selected off periods the unit can be restarted by pressing PROCESS. A few minutes operation should be allowed before use to allow water purity to reach a high.


Step 16 - Display viewing angle adjustment (screen 032)

The angle of the display can be adjusted up and down for better viewing of the screen.

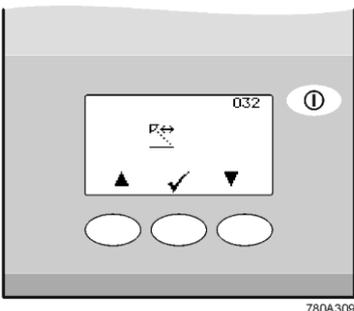

Viewing angle adjustment screen

1. PRESS and HOLD or to adjust the viewing angle
OR
PRESS to proceed to step 17.
2. PRESS .

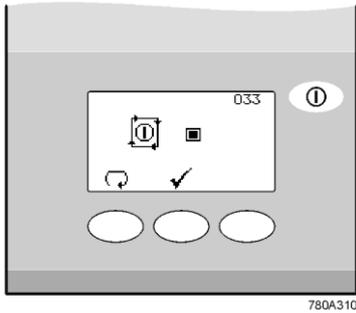


Auto-restart screen

Step 17 - Auto-restart (screen 033)

This allows selection of the AUTO restart option. If auto restart is selected the unit will automatically restart after a power failure. In manual mode the unit will stay in standby after a power failure.

1. PRESS to change mode (■ = ON)
OR
PRESS to proceed to step 18.
2. PRESS .

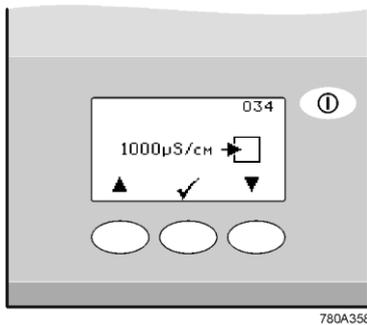


Feed Water Quality Screen

Step 18 - Feed Water Quality (screen 034)

An indication of the RO performance can be obtained using a calculation of ionic rejection in which the conductivity of the permeate is compared to that of the feedwater.

Enter the feedwater conductivity to obtain an accurate measurement.

1. PRESS to reduce the value to the correct reading
2. PRESS to accept

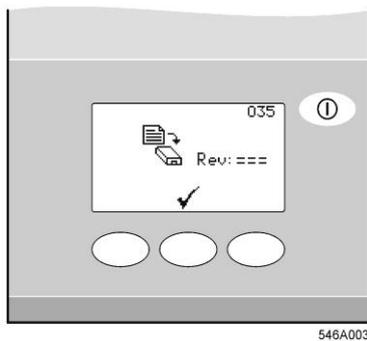


Data output screen

Step 19 - Data output (screen 035)

The **MEDICA EDI - US** can transmit operational data to an Xport storage device.

1. PRESS to format the storage device.
OR
PRESS to proceed to step 20.
2. PRESS .

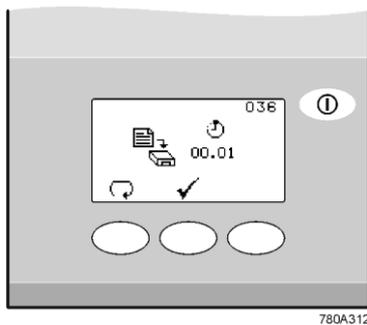

Data transmit screen

Step 20 - Data transmit (screen 036)


Setting the frequency of data transmittals to Xport.

1. PRESS to change transmit intervals (1, 5, 15, 30min / 1,6 hour)
OR
PRESS to proceed to step 21.
2. PRESS .

Note: PRESS PRINT during normal operation and current data will be transmitted if data output (035) is selected.

Step 21 - Programming of User Keys (screen 037)

During the life of the **MEDICA EDI - US** you may need to delete or add User keys to prevent or allow access to user choices. This feature is only available to the Master key holder.

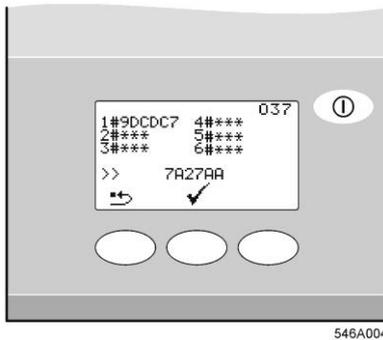

Programming of user keys screen




CAUTION! Do not press reset unless all User Keys are present.

1. PRESS ↵ to delete all User keys
2. PRESENT the new User key to the reader.
3. PRESS →** to load new User key identification.
4. REPEAT instruction 2 and 3 until all User keys are registered (maximum of 6 users).
5. PRESS ✓ to complete settings.

Step 22 - Programming of Sanitization Keys (screen 038)

Sanitization Keys are required to start a Sanitization. These keys can be deleted or added.

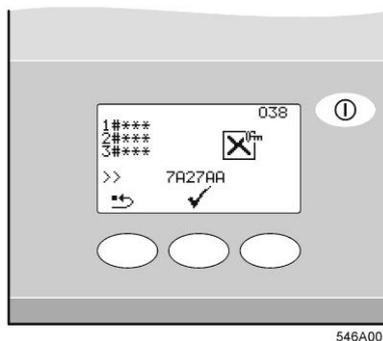

Programming of Sanitization Keys screen

1. PRESENT Master key to the reader
2. PRESS ✓ until the Sanitization Key programming screen appears.
3. PRESS ↵ to delete all Sanitization keys
4. PRESENT the new Sanitization Key to the reader.
5. PRESS →** to load new Sanitization key identification.
6. REPEAT instruction 4 and 5 until all Sanitization keys are registered (maximum of 2 users).
7. PRESS ✓ to complete settings.

3.3 Setting Consumable Replacement Reminders


Step 1 - Enter consumable replacement timer set-up

1. TURN unit off at electrical supply.

Note: Allow time for the display to go blank after the power is removed.

2. SWITCH electrical supply on.

3. PRESS Master Key and PRESS  button to enter the Consumable Timer set up.

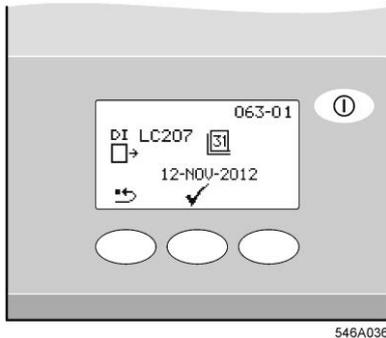

546A036


DI LC207 replacement date screen

Step 2 – Purification Pack (LC207) replacement date

1. PRESS  button to reset cartridge replacement date
OR
 PRESS  to accept replacement date and proceed to CVF Reminder.
2. PRESS  to confirm that resetting is required
OR
 PRESS *x* to abort reset.
3. PRESS .

Step 3 – Composite Vent filter (CVF) LC136M2 replacement date

1. PRESS  button to reset CVF replacement date
OR
 PRESS  to accept date and proceed to UV lamp Reminder.
2. PRESS  to confirm that resetting is required
OR
 PRESS *x* to abort reset.
3. PRESS .

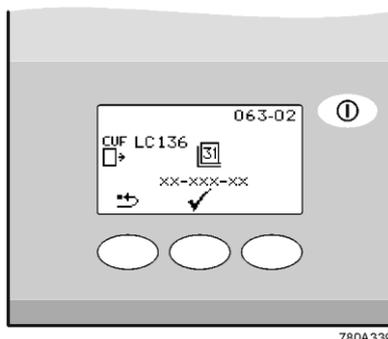


780A339

CVF replacement date screen


Step 4 - UV lamp (LC105) replacement date

1. PRESS  button to reset UV lamp replacement date
OR
 PRESS  to accept replacement date and proceed to Protek Filter Replacement.
2. PRESS  confirm that resetting is required
OR
 PRESS *x* to abort reset.
3. PRESS .

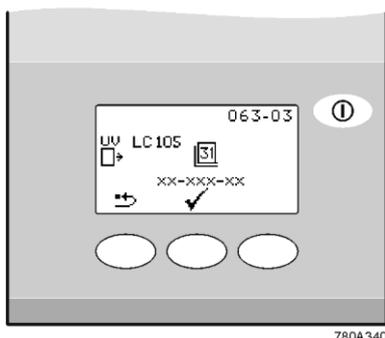


780A340


UV replacement date screen

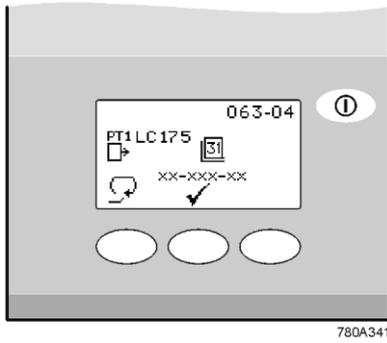


Pre-treatment cartridge replacement screen

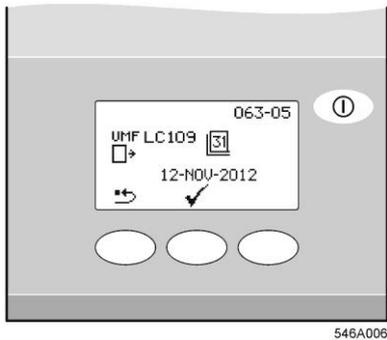

UMF replacement date screen

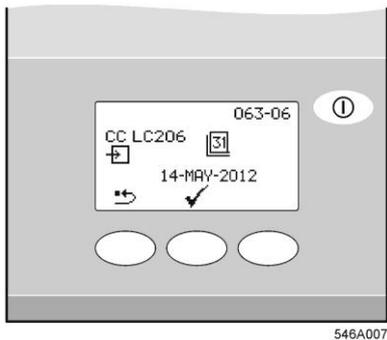

Conditioning Cartridge reminder screen

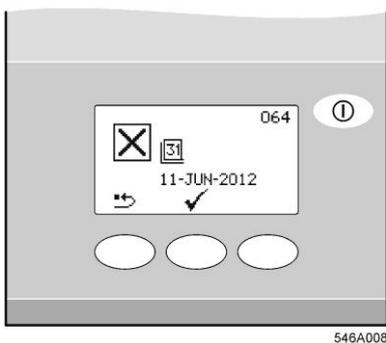

Sanitization reminder screen

Step 5 – Protek L2 (LC177).

1. PRESS ↵ button to reset Protek L2 (LC177) replacement date
OR
PRESS ✓ to accept replacement date and proceed to UMF replacement.
2. PRESS ✓ confirm that resetting is required
OR
PRESS X to abort reset.
3. PRESS ✓.

Step 6 – UMF (LC109) replacement date

1. PRESS ↵ button to reset UMF replacement date
OR
PRESS ✓ to accept replacement date and proceed to Sanitization reminder.
2. PRESS ✓ confirm that resetting is required
OR
PRESS X to abort reset.
3. PRESS ✓.

Step 8 – Conditioning Cartridge (LC206) reminder

1. Press ↵ button to reset cartridge replacement date
OR
Press ✓ to accept replacement date and complete settings.
2. PRESS ✓ confirm that resetting is required
OR
PRESS X to abort reset
3. PRESS ✓.

Step 7 - Sanitization reminder

1. PRESS ↵ button to reset Sanitization replacement date
OR
PRESS ✓ to accept replacement date and complete settings.
2. PRESS ✓ confirm that resetting is required
OR
PRESS X to abort reset.
3. PRESS ✓.

4. OPERATION

4.1 Night service/Operational day (screen 030)


Refer to Section 3.2 - Step 14 and 15.

Your **MEDICA EDI - US** can be programmed to operate on specific days between selected times. This is to optimize the efficiency of the unit and to minimize rises in water temperature.

During the 'sleep' period the unit will display .

It is possible to override this mode by PRESSING .

During the 'sleep' period the unit will run in intermittent recirculation (10 minutes every two hours) to maintain water purity around the distribution loop.

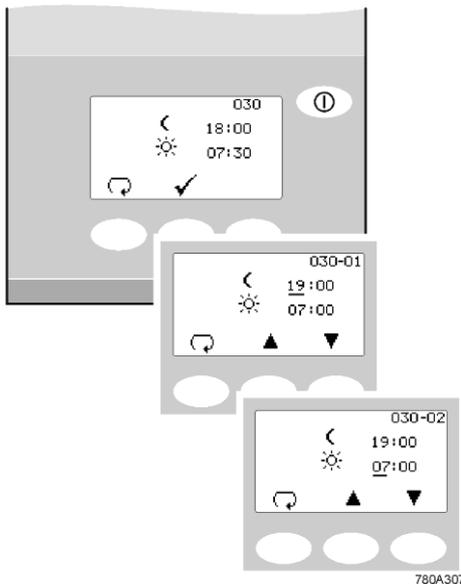

Night service/Operational day screen

4.2 Continuous Recirculation (24/7) (screen 029)

Refer to Section 3.2 - Step 13.

If the unit is set to continuous recirculation, it will constantly recirculate the water and fill the reservoir as required.

It is recommended that the system only runs in continuous mode when the demand for water is high (greater than 50% of the make up flow)

During recirculation the water temperature will tend to increase slowly.


5. MAINTENANCE

An approved supplier or distributor should carry out any maintenance work not included in this handbook.

Note: Disposal of all end of life consumable items should be in accordance with local statutory regulations



WARNING! ALWAYS CHECK THAT THE MAINS ELECTRICAL POWER AND FEED WATER ARE SWITCHED OFF BEFORE ATTEMPTING ANY MAINTENANCE PROCEDURE.

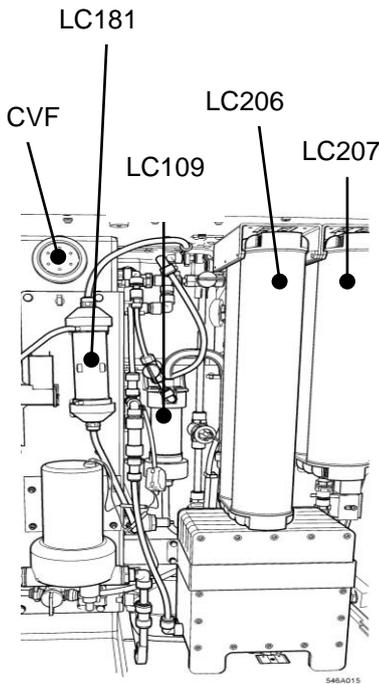

General View

5.1 Replacing Composite Vent Filter (L136)

The Composite Vent Filter (CVF) should be replaced in the following circumstances:

- When indicated by the alarm (screen prompt) or after a maximum of six months.
1. ENSURE process is OFF and isolate power.
 2. OPEN front doors and locate CVF.
 3. UNSCREW old CVF and discard
 4. UNPACK new CVF.
 5. WRITE the installation date on the label of the filter for future reference.
 6. INSTALL filter.
 7. SWITCH on power.
 8. RESET consumable reminder as described in Section 3.3, Step 3.

5.2 Replacing Purification Cartridge LC207 or Conditioning Cartridge LC206.

• LC206:

- When the purity of water from the unit starts to deteriorate.


When indicated by the consumable alarm or after a maximum of six months.

1. ENSURE process is OFF and isolate power.
2. OPEN right hand door and LOCATE LC207 cartridge, the shorter pack on the right (above the EDI pulse module)
3. PUSH Purification cartridge FORWARD, then LIFT, finally PULL to remove the used cartridge. Discard.
4. UNPACK new LC207.
5. REMOVE sealing plugs from inlet and outlet ports.
6. Wet 'O' rings and SLIDE new cartridge into position pushing upwards against pack reader contacts.
7. EASE back and ENSURE that the pack is fully engaged (down) in the retainers.
8. SWITCH on power. FOLLOW screen prompts to accept new cartridge replacement date
9. START unit and allow to circulate until water quality is achieved.
10. RESET reminder as described in section 3.3, Step 2.

• LC207:

- The replacement of the LC206 pack is dictated by water usage and RO water quality.
- Pack needs replacing when indicated by alarm 57.
- Alarm 56 will come up before, when 10% of pack life remains, and will indicate approximate remaining life in days, based on recent usage patterns

1. ENSURE process is OFF and isolate power.
2. OPEN right hand door and LOCATE LC206 cartridge, the longer pack on the left (above the EDI module).
3. PUSH purification cartridge FORWARD, then LIFT, finally PULL to remove the used cartridge. Discard.
4. UNPACK new LC206.
5. REMOVE sealing plugs from inlet and outlet ports.
6. Wet 'O' rings and SLIDE new cartridge into position pushing upwards against pack reader contacts.
7. EASE back and ENSURE that the pack is fully engaged (down) in the retainers.
8. SWITCH on power. FOLLOW screen prompts to accept new cartridge replacement date.
9. RESET reminder as described in section 3.3, Step 2.

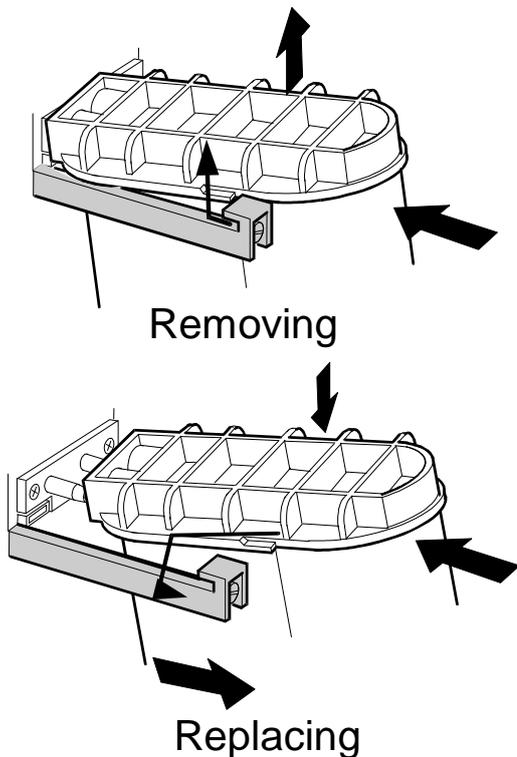

Purification Cartridge (LC207) and Conditioning Cartridge (LC206) removal and replacement.

5.3 Replacing Ultra-Microfilter (LC109)


The Ultra-Microfilter (UMF) should be replaced in the following circumstances:

- When indicated by the consumable alarm
- After a maximum of 6 months use
- Product flow is no longer adequate
- Product water bacterial content is too high.

1. ISOLATE power
2. OPEN front doors and LOCATE UMF
3. UNSCREW lower connection


Note: At this stage a small quantity of water may spill.

4. UNSCREW UMF from top connection
5. INSTALL new UMF
6. RECONNECT
7. RINSE water to drain to recommission the system.
8. RESET reminder as described in section 3.3, Step 6.

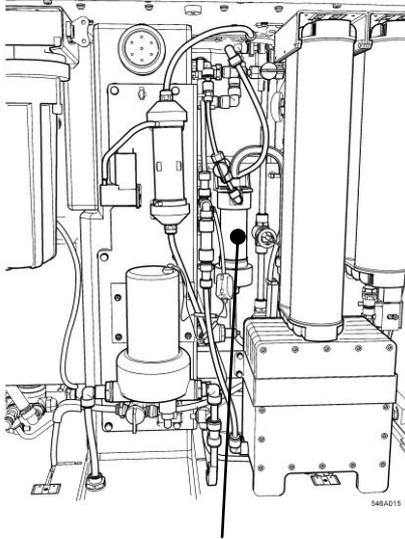

Ultra-Microfilter

Ultra-Microfilter position

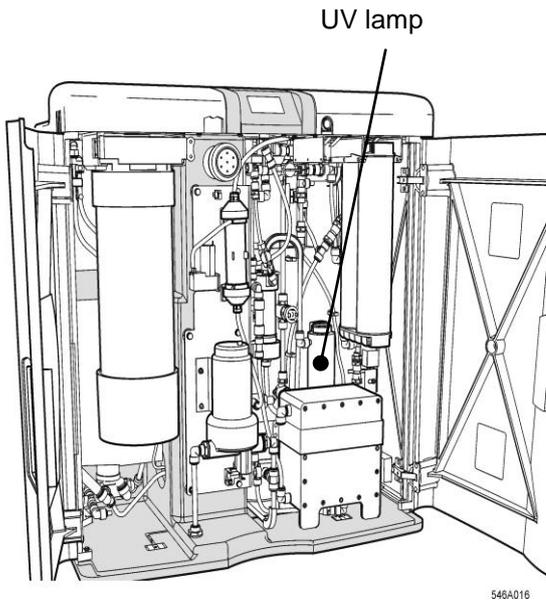

UV lamp

UV lamp position

5.4 Replacing UV lamp (LC105)

The UV lamp should be replaced in the following circumstances:

- When indicated by the consumable alarm.
 - After a maximum of twelve months use.
1. ENSURE process is OFF and ISOLATE power.
 2. OPEN the front doors.
 3. REMOVE pack to improve access to the UV lamp assembly
 4. LOCATE UV assembly on right-hand side of unit.
 5. REMOVE retaining spring clips from top and bottom of the lamp assembly
 6. REMOVE electrical connectors from top and bottom and take lamp from housing.




WARNING! IT IS STRONGLY RECOMMENDED THAT DURING THE HANDLING OF THE LAMP CUT-RESISTANT GLOVES BE WORN.

7. DISCARD lamp.




Caution! Dispose of lamp in accordance with local authority regulations.

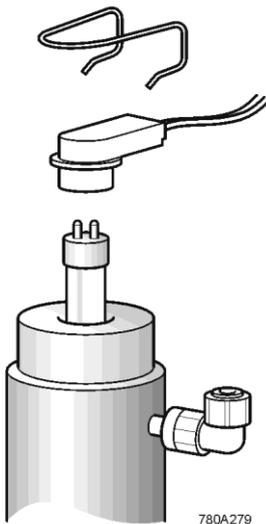


UV Lamp connection

8. REMOVE new lamp from packaging and follow instructions included for cleaning.
9. REFIT into UV assembly.
10. RECONNECT to electrical connectors top and bottom.
11. REFIT retaining spring clips top and bottom
12. REPLACE Medpure pack (if fitted).
13. RESET consumable reminder as described in Section 3.3, Step 4.


Protek L2

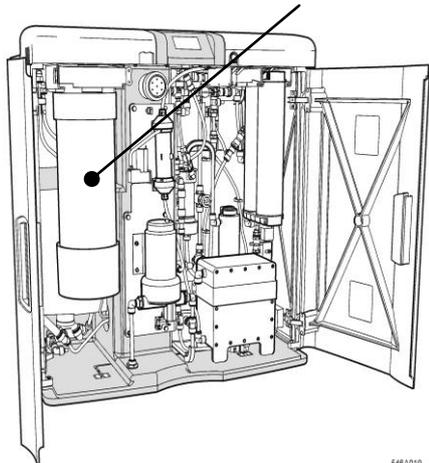

Protek position

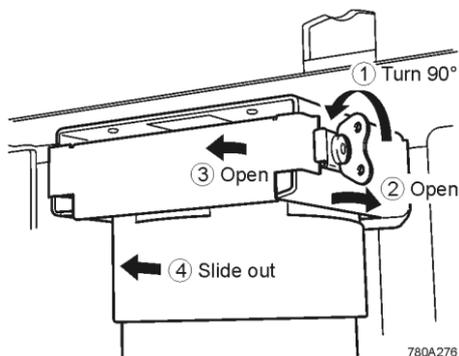

Protek L2 (LC177)

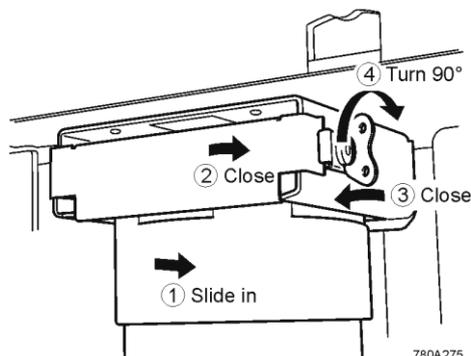

Protek L2 (LC177)


5.5 Replacing Protek L2 (LC177) or Carbon Filter (LC123) Pre-treatment.

The replacement frequency of the pre-treatment cartridge is dictated by the purity of the feed water. It should be replaced in the following circumstances:

- When indicated by the consumable alarm
- When indicated by the alarm 75
- After replacement of RO modules

5.5.1 Protek L2 (LC177)

1. ENSURE process is OFF and ISOLATE power.
2. OPEN front doors
3. LOCATE Protek
4. Remove lid from break tank and RELEASE any residual system pressure by OPERATING float valve.
5. TWIST to RELEASE clip and OPEN retaining clamp at top of cartridge Protek L2 (LC177).
6. PULL Protek L2 (LC177) forward
7. If replacing a Protek, discard the used cartridge in accordance with local regulations.
8. REMOVE new cartridge from packaging.
9. REMOVE sealing plugs from inlet and outlet ports.
10. Wet 'O' rings and SLID Protek until it is fully engaged
11. CLOSE retaining clamp and fasten CLIP
12. The unit will automatically recognize the new cartridge and date.
13. RESET reminder as described in section 3.3, Step 5.

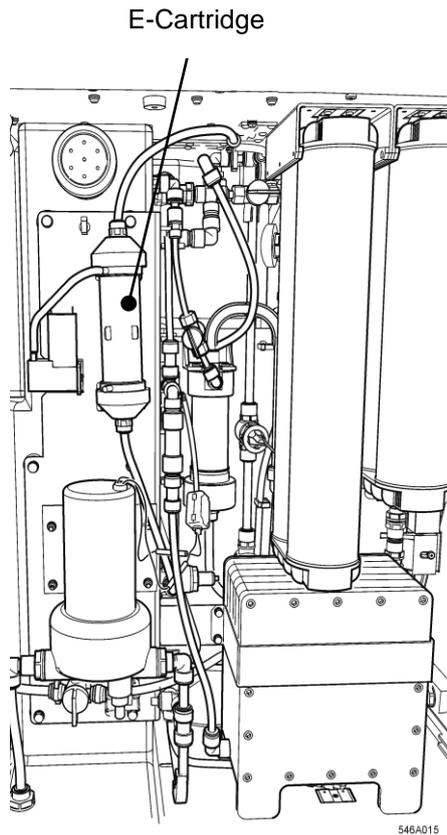

5.6 Replacing E- Cartridge (LC181)

The E- cartridge should be replaced:

- When indicated by the alarm
1. ENSURE process is OFF and ISOLATE power.
 2. OPEN front doors
 3. LOCATE E – cartridge
 4. DISCONNECT tubing to air pump and inlet and outlet tubing from the E cartridge.
 5. REMOVE used cartridge from clip and discard in accordance with local regulations.
 6. TAKE new cartridge
 7. REMOVE sealing plugs from inlet and outlet ports
 8. LOCATE new cartridge in clip
 9. RE-CONNECT tubing

5.7 RO modules

The reverse osmosis modules should be replaced if permeate water purity or flow rate is not adequate and does not meet predicted or previous performance. For information regarding the replacement of the reverse osmosis module contact Customer Services.

6. SANITIZATION

6.1 Liquid sanitization

The RO is sanitized to reduce the growth of microbiological contamination within the RO module. The **MEDICA** has a built-in reservoir which allows the sterilant to be introduced into the water feeding the RO. Please read this entire section to become familiar with the procedure before you start.

The sanitization frequency required is dependent on the feedwater, local environment, usage patterns and application. As a general rule, ELGA LabWater recommends sanitization of the RO at 6 monthly intervals. However, the period between sanitizations could be extended in particular circumstances. For example, microbial growth will usually be lower with clean feed water, well-maintained pre-treatment, low temperatures and heavy usage.



WARNING! ALWAYS WEAR RUBBER GLOVES, APRON AND FACEMASK. DO NOT BREATHE FUMES OR ALLOW TO COME IN CONTACT WITH SKIN OR EYES. ALWAYS FOLLOW RECOMMENDATIONS FOUND IN THE MANUFACTURERS MATERIAL SAFETY DATA SHEET AND ANY APPLICABLE OSHA STANDARDS FOR THE CHEMICAL BEING USED.




WARNING! LABEL THE MACHINE WITH APPROPRIATE WARNING SIGNS SUCH AS “DO NOT USE/ CONTAINS STERILANT” (NOT PROVIDED)

6.1.1 RO sanitization (Step 1) RO sanitization.


Sanitization of the finished installation is essential to ensure that the system is properly commissioned and capable of achieving the bacteriological control required for a typical application.

1. PRESENT GREEN sanitization key to reader.
2. PRESS ✓ to enter sanitization.

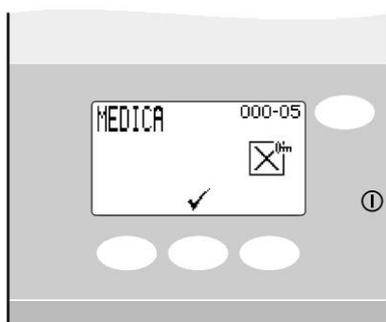


780A449

Sanitization selection screen

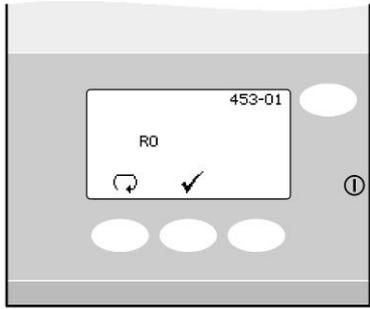

780A430

Enter RO sanitization

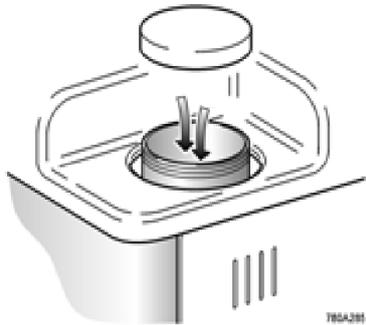


780A431

Large white cap removal to insert sanitant

3. PRESS ✓ to enter RO sanitization.

- 4. UNSCREW the **large** white lid on the left of the unit.
- 5. Measure 50ml of Minncare cold sterilant and slowly pour into the break tank.
- 6. REPLACE lid.

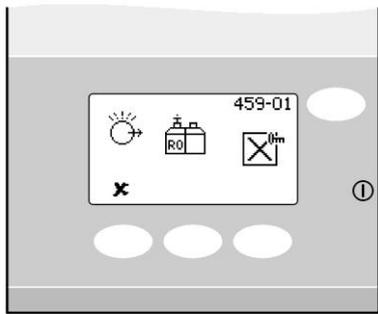


780A433

Confirmation chemical is added

7. PRESENT GREEN sanitization to reader to confirm that the chemical has been added.


- 8. The system now enters a period of automatic sanitization lasting approximately 45 minutes

Once the step within the sanitization is complete the unit will go straight into process and continue to fill the recirculation reservoir

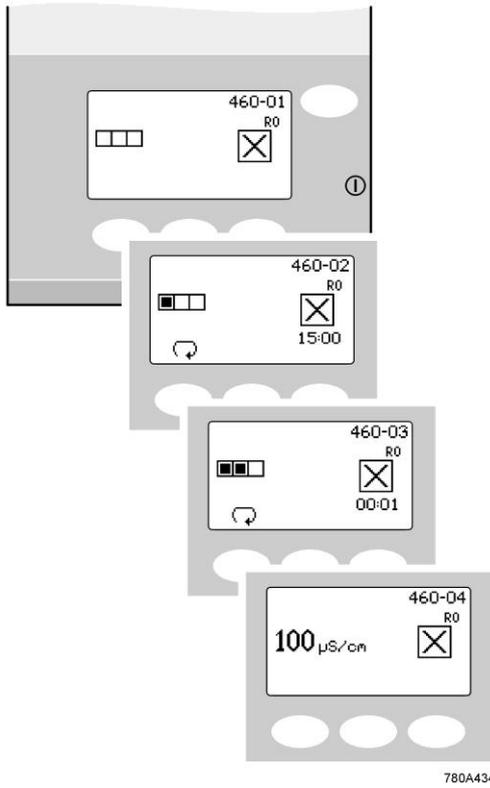

RO sanitization process steps

- 9. PRESS ✓ to accept the sanitization reminder

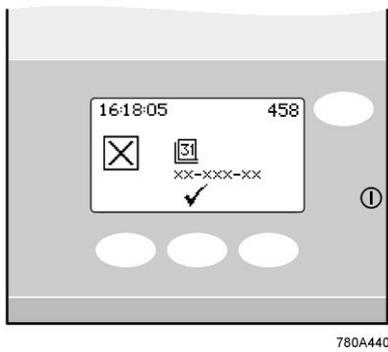

Set Sanitization reminder

6.2 Tablet sanitization

The RO is sanitized to reduce the growth of microbiological contamination within the RO module. The **MEDICA** has a built-in reservoir which allows the sterilant to be introduced into the water feeding the RO. Please read this entire section to become familiar with the procedure before you start.

The sanitization frequency required is dependent on the feedwater, local environment, usage patterns and application. As a general rule, ELGA LabWater recommends sanitization of the RO at 6 monthly intervals. However, the period between sanitizations could be extended in particular circumstances. For example, microbial growth will usually be lower with clean feed water, well-maintained pre-treatment, low temperatures and heavy usage.



WARNING! ALWAYS WEAR RUBBER GLOVES, APRON AND FACEMASK. DO NOT BREATHE FUMES OR ALLOW TO COME IN CONTACT WITH SKIN OR EYES. ALWAYS FOLLOW RECOMMENDATIONS FOUND IN THE MANUFACTURERS MATERIAL SAFETY DATA SHEET AND ANY APPLICABLE OSHA STANDARDS FOR THE CHEMICAL BEING USED.




WARNING! LABEL THE MACHINE WITH APPROPRIATE WARNING SIGNS SUCH AS “DO NOT USE/ CONTAINS STERILANT” (NOT PROVIDED)

6.2.1 RO sanitization (Step 1).

Sanitization of the finished installation is essential to ensure that the system is properly commissioned and capable of achieving the bacteriological control required for a typical application.


1. PRESENT GREEN sanitization key to reader.
2. PRESS ✓ to enter sanitization.

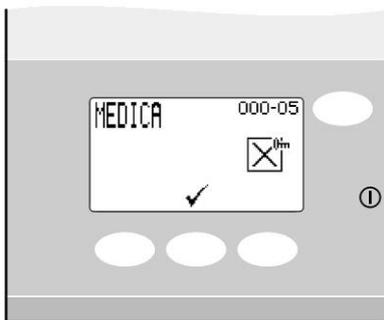

780A449

Sanitization selection screen

- 3. PRESS ✓ to enter RO sanitization.

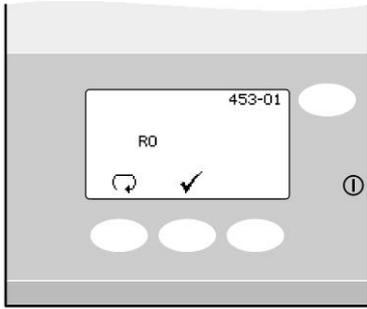


780A430

Enter RO sanitization




CAUTION! ENSURE THE CORRECT DOSAGE OF CHEMICAL IS USED.

- 4. UNSCREW the **large** white lid on the left of the unit.
- 5. CUT one EfferSan™ tablet into four quarters using tablet cutter.

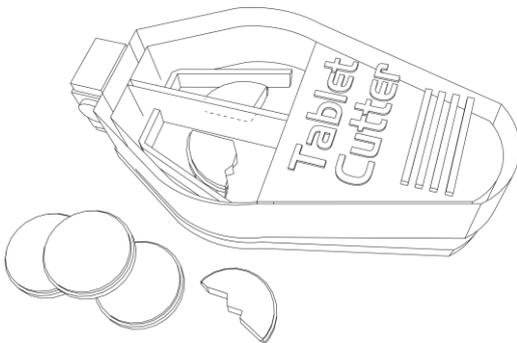


Tablet cutter for cutting EfferSan™ tablets.

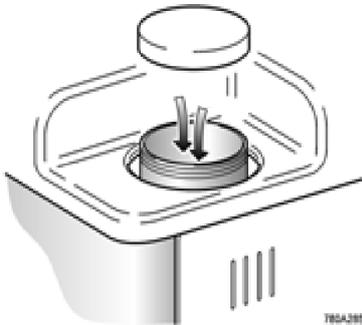

780A430

Large white lid removal to insert sanitant.

- 6. INSERT ¼ of EfferSan™ tablet into the break tank.
- 7. REPLACE lid.

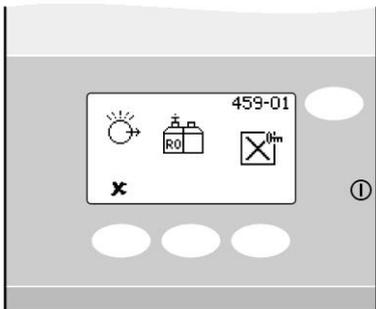


780A433

Confirmation chemical is added

- 8. PRESENT GREEN sanitization to reader to confirm that the chemical has been added.

- 9. The system now enters a period of automatic sanitization lasting approximately 45 minutes


Once the step within the sanitization is complete the unit will go straight into process and continue to fill the recirculation reservoir

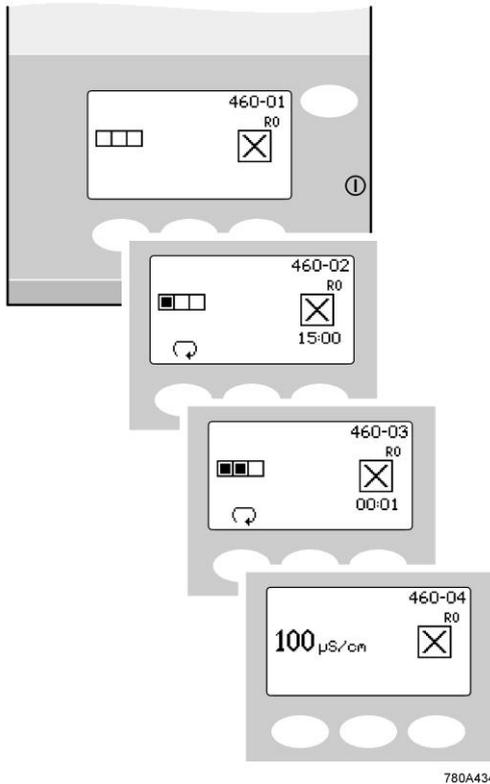

780A434

RO sanitization process steps

- 10. PRESS ✓ to accept the sanitization reminder

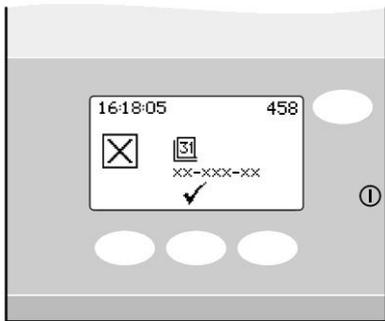


780A440


Set Sanitization reminder

7. EMERGENCY BYPASS

Contact your local service provider before putting the unit into emergency bypass.

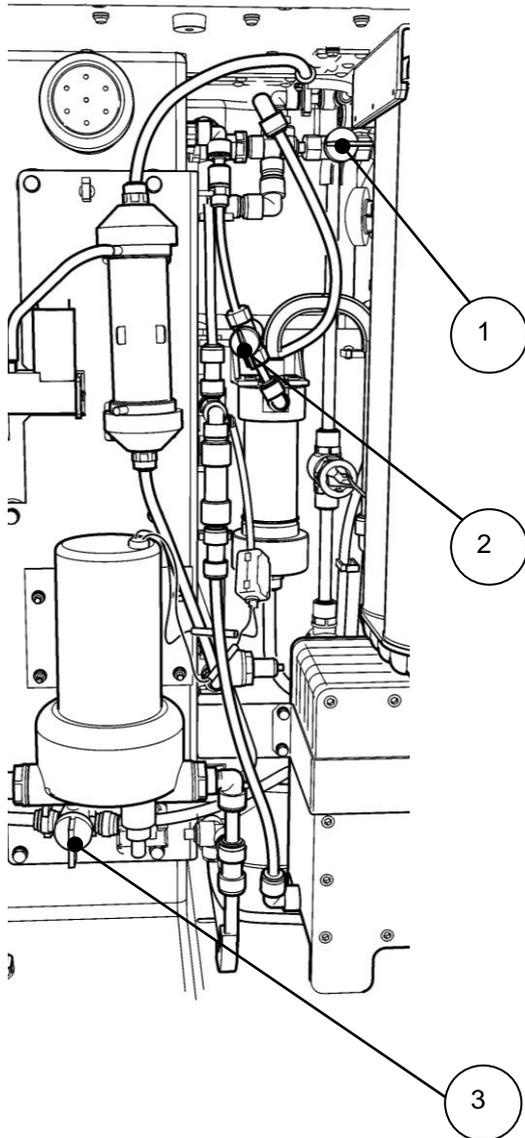

7.1 Emergency bypass

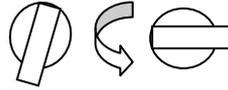1. Ensure the process is OFF
2. Disconnect the electrical supply
3. **Open the doors.**
4. Close valve 1 
5. Close valve 2 
6. **Remove yellow plastic cover 3**
7. Open valve 3 
8. Close valve 3 when there is no requirement for purified water.
9. **Monitor water purity independently or install new LC207 cartridges to guarantee water purity.**

7.2 Calculation for DI Cartridge (LC207) Cartridge life

Calculate the maximum life of LC207 under emergency bypass processing and with your mains water conditions to maintain a supply purity of 1MΩ.cm to the analyzer

The maximum life expectancy of the LC207 cartridge in emergency bypass mode while maintaining a supply purity of 1MΩ at the analyzer, may be calculated as follows:

Single Pack

$$\text{Time between pack changes} = \frac{60,000}{\text{Conductivity of mains water} \times \text{Water consumption of analyzer (l/hr)}}$$

Example

Mains water conductivity = 605 μS/cm

Water consumption of analyzer = 20 l/hr

$$\text{Time between pack changes} = \frac{120,000}{605 \times 20}$$

Time between pack changes = 4.9 hr

In the absence of information on mains water conductivity, assume a pack life of 2 hr if the analyzer consumes 20 l/hr. If flows start to decline this may be improved by changing the UMF.

The capacity will be less (typically 80% less) if a water purity of 10 MΩ.cm is needed.

8. TROUBLE SHOOTING

If a problem occurs, the unit will normally sound an alarm and the respective icon will flash. The audible alarm can be silenced by pressing the mute button.

The appropriate action to then take (for the more common alarm displays) is shown on the Quick Reference Guide (fitted to the inside of the door of the **MEDICA EDI - US.**) If the unit cannot be repaired using the information shown, please call your local ELGA LabWater representative (See Section 12.0 - Useful Contact Details).



WARNING! ALWAYS ENSURE THAT THE MAINS POWER SUPPLY IS ISOLATED BEFORE WORKING INSIDE THE UNIT.

Alarm Summary

Display id	icon	Alarm	Recommended Action
672-98		Leak Detection	<ul style="list-style-type: none"> PRESS mute to silence alarm REMOVE power rectify leak and dry sensors Dry the contacts POWER on to reset.
672-97		Reservoir level controls	<ul style="list-style-type: none"> PRESS mute to silence alarm REMOVE power Contact Technical Support
672-96		High loop pressure	<ul style="list-style-type: none"> PRESS mute to silence alarm REMOVE power Contact Technical Support
672-95		Permeate over temperature	<ul style="list-style-type: none"> PRESS mute to silence alarm Potable feed temperature is unacceptably high >50°C Contact Technical support
672-94		Recirc. loop over temperature	<ul style="list-style-type: none"> PRESS mute to silence alarm Water temperature is unacceptably high >50°C Contact Technical support
-93		Protek L2 incorrectly installed	<ul style="list-style-type: none"> PRESS mute to silence alarm CONFIRM Protek is correctly installed POWER OFF/ ON to reset
-92		LC207 incorrectly installed	<ul style="list-style-type: none"> PRESS mute to silence alarm CONFIRM Medpure L1 is correctly installed POWER OFF/ ON to reset
-90		Water purity alarm	<ul style="list-style-type: none"> PRESS mute to silence alarm CONFIRM purity displayed is suitable for application. CHANGE LC207 if water purity is insufficient
-89		Water temperature alarm	<ul style="list-style-type: none"> PRESS mute to silence alarm CONFIRM temperature is suitable for application. ADJUST alarm setpoint or DIRECT water to drain to Introduce cool water
-88		Permeate purity alarm	<ul style="list-style-type: none"> PRESS mute to silence alarm CONFIRM purity displayed is suitable for application. ADJUST alarm setpoint or contact local service provider.
-87		Permeate temperature alarm	<ul style="list-style-type: none"> PRESS mute to silence alarm CONFIRM temperature is suitable for application. ADJUST alarm setpoint or confirm potable supply temperature is suitable
-86		Break tank low	<ul style="list-style-type: none"> PRESS mute to silence alarm. CONFIRM pressure and flow of potable supply If problem persists contact local service provider
-85		UV lamp fail	<ul style="list-style-type: none"> PRESS mute to silence alarm CHANGE lamp at the next opportunity
60		Invalid C.C.	<ul style="list-style-type: none"> Check if you have re-installed old LC206 Put in new LC206
57		Replace LC206	<ul style="list-style-type: none"> LC206 is expired. Replace with new one.
56		LC206 reminder	<ul style="list-style-type: none"> Make sure you have a new LC206 ready in stock for replacement 10% of life left.

9. KEY TO CONTROL PANEL

ICON	DESCRIPTION	ICON	DESCRIPTION
	Accept		User Recognition Key
	Auto restart		Locked
	Scroll back		Menu
	Step back		Mute alarm
	Bell		Night
	Standby		Option OFF
	Calibration point		Option ON
	Cancel		Output
	Clock		Recirculate
	Date		Reset
	Day		Right
	Down		Sanitization PASSkey
	Drain		Sanitization reminder
	Save data		Scroll
	Hazard		Transport mode
	Feed		UP
	Fill		Viewing angle
	Add chemical		Recirc

Dimensions		
Height	820 mm (32.8")	834 mm (33.4") including castors
Width	794 mm (31.8")	
Depth	470 mm (18.8")	
Unit Type	MEDICA EDI 50	MEDICA EDI 100
Supply weight	58 Kg	65 Kg
Operational weight	115 Kg	122 Kg
Installation	Floor/ bench	

Connections	
Inlet	10 mm OD tube
Drain	10 mm OD tube
EDI drain	10 mm OD tube
Recirculation loop outlet*	12 mm OD tube
Recirculation loop inlet*	12 mm OD tube

Electrical Requirements	
Mains input	230 V ac, 50 Hz 115 V ac, 60 Hz
System control voltage (not including pumps and UV)	24 V dc
Power consumption (peak demand)	650 VA
Electrical protection rating	10 amp

Safety Features
Power fail safe
Water temperature alarm*
Water purity alarm*
Leak detection alarm
Access restricted by User Recognition key
Low voltage control circuit - 24 V dc
Visual alarms
Audible alarms

* Output to operate 24 V dc relay.

Product Water Specification	
Flowrate	MEDICA EDI
Daily usage - typical	960
Daily usage - maximum	2880
Inorganic	>10 MΩ. cm @ 25°C
TOC ppb	<30
Bacteria	<1 CFU/ml**
Particles	0.05 µm

** System to be regularly sanitized and installed following ELGA LabWater installation design guidelines.

As part of our policy of continual improvement, we reserve the right to alter the specifications given in this document.

11. WARRANTY/CONDITIONS OF SALE

ELGA LabWater is a trading name of VWS (UK) Ltd.

General Limited Warranty

VWS (UK) Ltd. warrants the products manufactured by it against defects in materials and workmanship when used in accordance with applicable instructions for a period of one year from the date of shipment for the products. VWS (UK) Ltd. MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED. THERE IS NO WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. The warranty provided herein and the data, specifications and descriptions of the VWS (UK) Ltd. products appearing in VWS (UK) Ltd. published catalogues and product literature may not be altered except by express written agreement signed by an officer of VWS (UK) Ltd. Representations, oral or written, which are inconsistent with this warranty or such publications are not authorized and, if given, should not be relied upon.

In the event of a breach of the foregoing warranty, VWS (UK) Ltd. sole obligation shall be to repair or replace, at its option, any product or part thereof that proves to be defective in materials or workmanship within the warranty period, provided the customer notifies VWS (UK) Ltd. promptly of any such defect. The exclusive remedy provided herein shall not be deemed to have failed of its essential purpose so long as VWS (UK) Ltd. is willing and able to repair or replace any nonconforming VWS (UK) Ltd. product or part. VWS (UK) shall not be liable for consequential, incidental, special or any other indirect damages resulting from economic loss or property damage sustained by any customer from the use of its products.

VWS (UK) Ltd. Warranty

VWS (UK) Ltd. warrants the water systems manufactured by it, BUT EXCLUDING MEMBRANES AND PURIFICATION PACKS, against defects in materials and workmanship when used in accordance with the applicable instructions and within the operating conditions specified for the systems for a period of one year from the earlier of:

- a) the date of installation, or
- b) the 120th day following the date of shipment.

VWS (UK) LTD. MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED. THERE IS NO WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. The warranty provided herein and the data, specifications and descriptions of the VWS (UK) Ltd. systems appearing in VWS (UK) Ltd. published catalogues and product literature may not be altered except by express written agreement signed by an officer VWS (UK) Ltd. Representations, oral or written, which are inconsistent with this warranty or such publications are not authorised and, if given, should not be relied upon.

In the event of a breach of the foregoing warranty, VWS (UK) Ltd. sole obligation shall be to repair or replace, at its option, any product or part thereof that proves to be

defective in materials or workmanship within the warranty period, provided the customer notifies VWS (UK) Ltd. promptly of any such defect. The cost of labor for the first ninety (90) days of the above warranty period is included in the warranty; thereafter, labor cost shall be at the customer's expense. The exclusive remedy provided herein shall not be deemed to have failed of its essential purpose so long as VWS (UK) Ltd. is willing and able to repair or replace any nonconforming VWS (UK) Ltd. system or component part. VWS (UK) Ltd. shall not be liable for consequential, incidental, special or any other indirect damages resulting from economic loss or property damage sustained by any customer from the use of its process systems.

Products or components manufactured by companies other than VWS (UK) Ltd. or its affiliates ("Non- VWS (UK) Ltd. products") are covered by the warranty, if any, extended by the Product manufacturer. VWS (UK) Ltd. hereby assigns to the purchaser any such warranty; however VWS (UK) LTD. EXPRESSLY DISCLAIMS ANY WARRANTY WHETHER EXPRESSED OR IMPLIED, THAT THE NON - VWS (UK) Ltd. PRODUCTS ARE MERCHANTABLE OR FIT FOR A PARTICULAR PURPOSE.

NOTICE

VWS (UK) Ltd. is constantly striving to improve its products and services. Consequently, the information in this document is subject to change without notice and should not be construed as a commitment by VWS (UK) Ltd. Also VWS (UK) Ltd. assumes no responsibility for any errors that may appear in this document. This manual is believed to be complete and accurate at the time of publication. In no event shall VWS (UK) Ltd. be liable for incidental or consequential damages in connection with or arising from the use of this manual.

VWS (UK) Ltd. warrants its products against defects in materials and workmanship as described in the Warranty statement on the preceding pages.

12. USEFUL CONTACT DETAILS

ELGA LabWater

Lane End Industrial Park

High Wycombe

Bucks HP14 3BY

UK

Tel: +44 (0) 203 567 7300

Fax: +44 (0) 203 567 7205

E-mail: techsupport@elgalabwater.com

For the address of the nearest ELGA LabWater Sales and Service office visit the country list on our website.

<http://www.elgalabwater.com>

or contact ELGA at the number above.

ELGA® is the global laboratory water brand name of Veolia Water.

ELGA®, PURELAB®, MEDICA® and Centra® are trademarks of ELGA LabWater (VWS UK Ltd).

Minnicare is a trademark of Minntech Corporation

©2007-2011 ELGA LabWater/VWS (UK) Ltd. All rights reserved

Copyright Note

The information contained in this document is the property of VWS (UK) Ltd. and is supplied without liability for errors or omissions.

No part of this document may be reproduced or used except as authorized by contract or other written permission from VWS (UK) Ltd. The copyright and all restrictions on reproduction and use apply to all media in which this information may be placed.

VWS (UK) Ltd. pursue a policy of continual product improvement and reserve the right to alter without notice the specification, design, price or conditions of supply of any product or service.