PURELAB[®]

Chorus 2

Operator Manual
MANU40003 VER 03 3/19

WATER PURIFICATION



ELGA  VEOLIA

WELCOME

Thank you for purchasing the PURELAB® Chorus 2. Please read the following Health and Safety notice before proceeding.

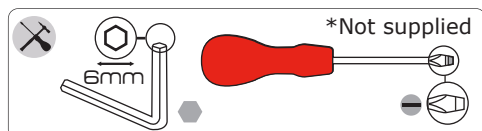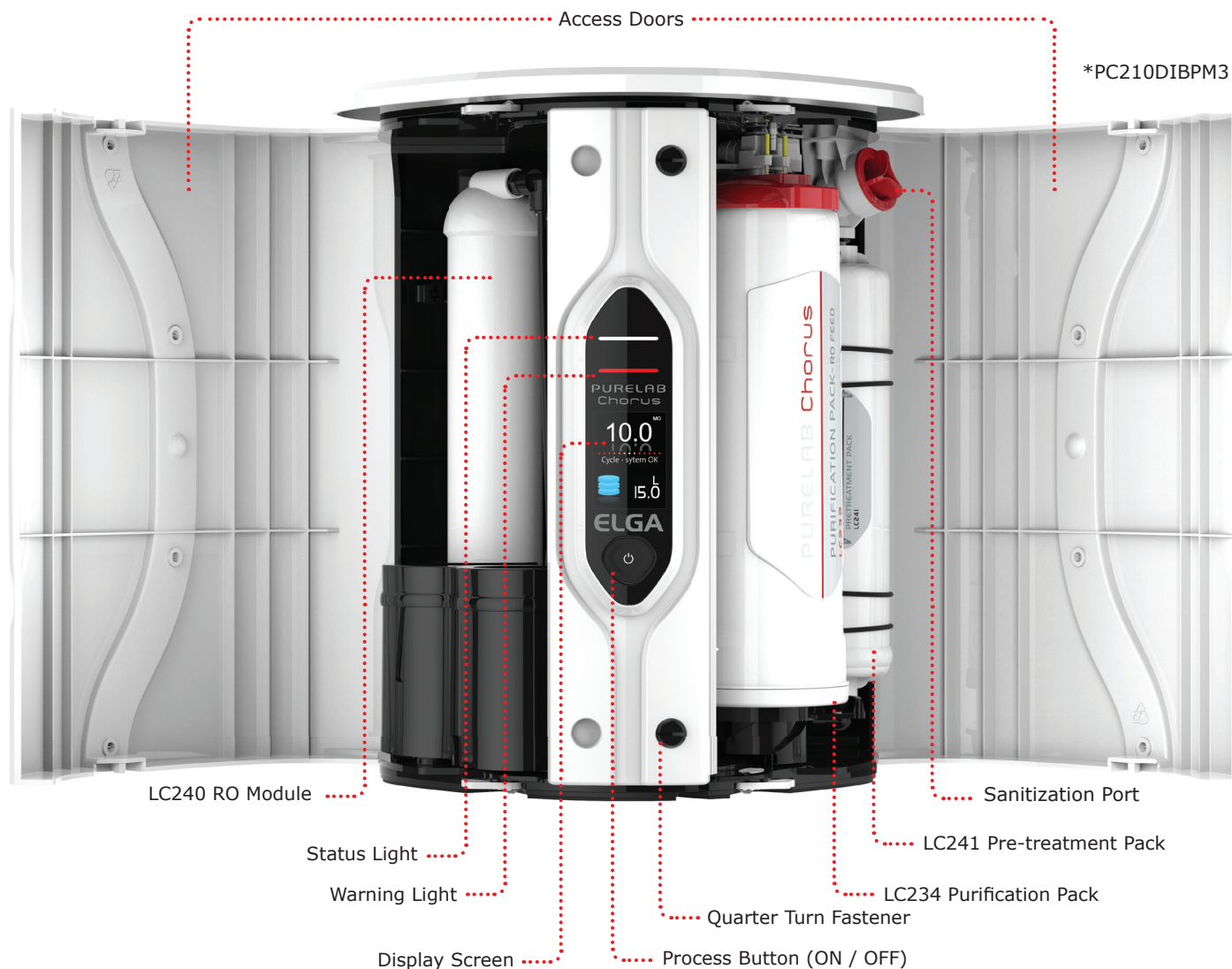Use a slotted screw driver to release quarter turn fastener slots on Chorus.

HEALTH AND SAFETY

Please follow the safety information detailed in this instruction sheet. Failure to observe these instructions could result in damage to the equipment, processes, associated equipment and injury or death to persons.

1: Electricity

WARNING!

POSITION THE POWER SUPPLY SO THAT IT CANNOT COME INTO CONTACT WITH WATER.

2: Residual Preservative

WARNING!

DURING THE COMMISSIONING CYCLE TRACE RESIDUAL PRESERVATIVE IS FLUSHED FROM THE SYSTEM. A MATERIAL SAFETY DATA SHEET IS AVAILABLE UPON REQUEST. techsupport@elgalabwater.com


Control of Substances Hazardous to Health (COSHH)

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

Example One: 1 x PURELAB® Chorus 2 (RO/DI) + 15L or 30L Reservoir LA757 / LA578

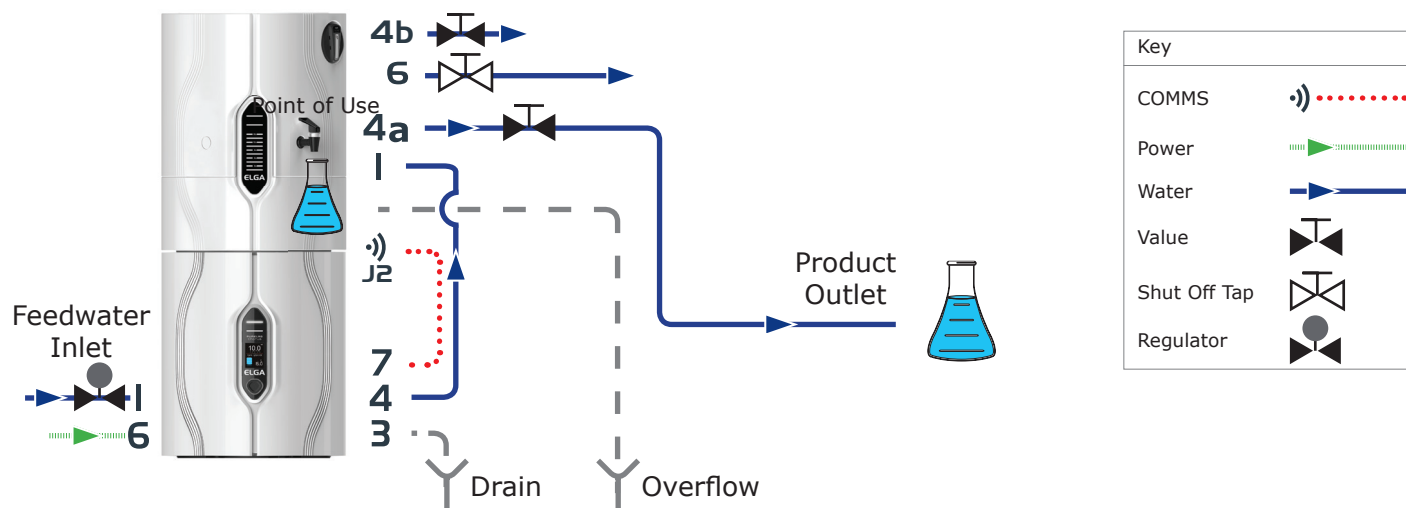Example Two: 2 x PURELAB® Chorus 2 (RO/DI) + 60L Reservoir LA759

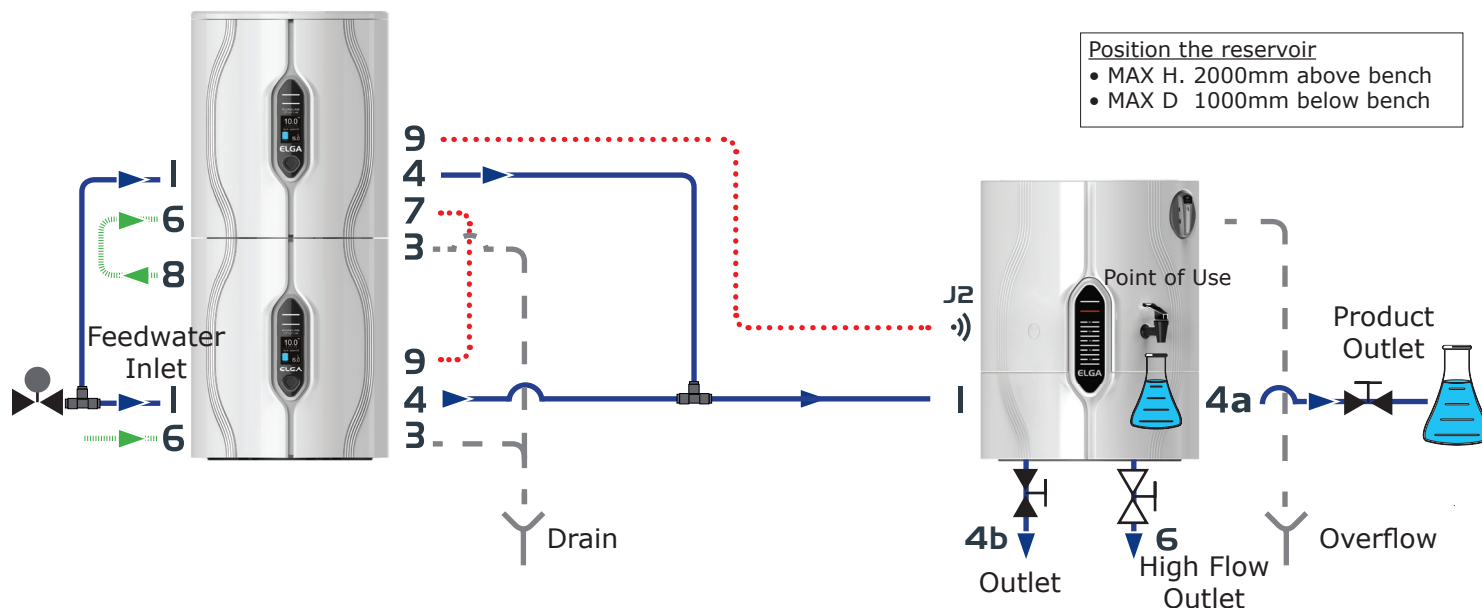
Example Three: 4 x PURELAB® Chorus 2 (RO/DI) + 100L Reservoir LA760

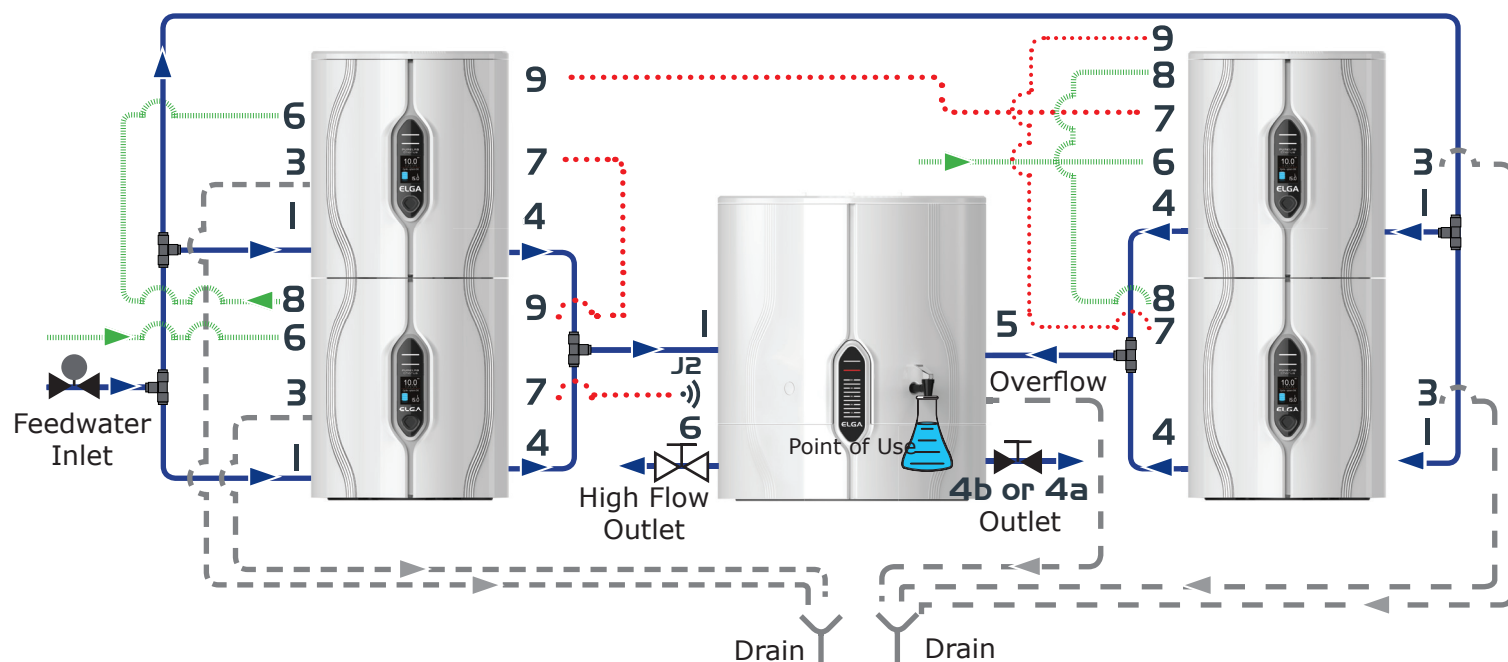SYSTEM CONFIGURATION

Example Four: 1 x PURELAB[®] Chorus 2 (RO/DI) + 15L or 30L Reservoir LA757 / LA578

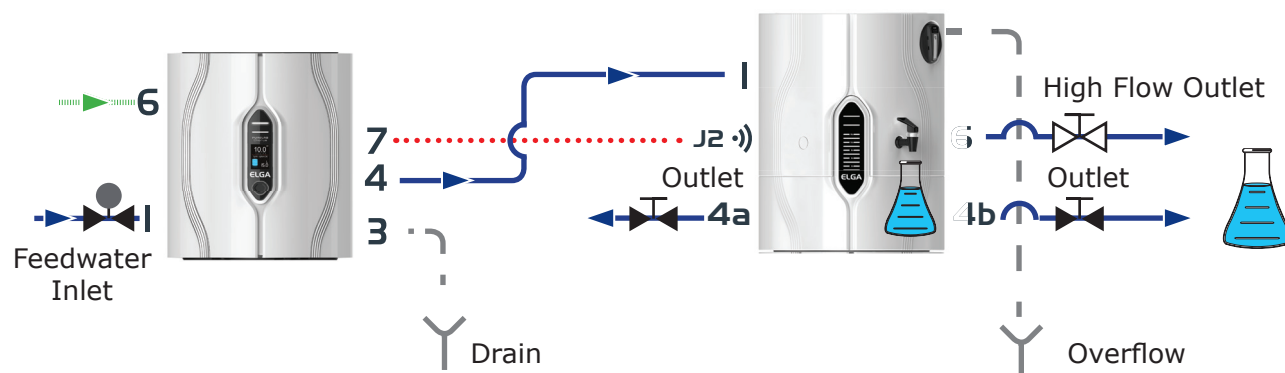
Further advice on configuration connections.

SET - UP

Ia MODULAR STACKING BENCH MOUNTED RESERVOIR ON TOP (OPTIONAL)

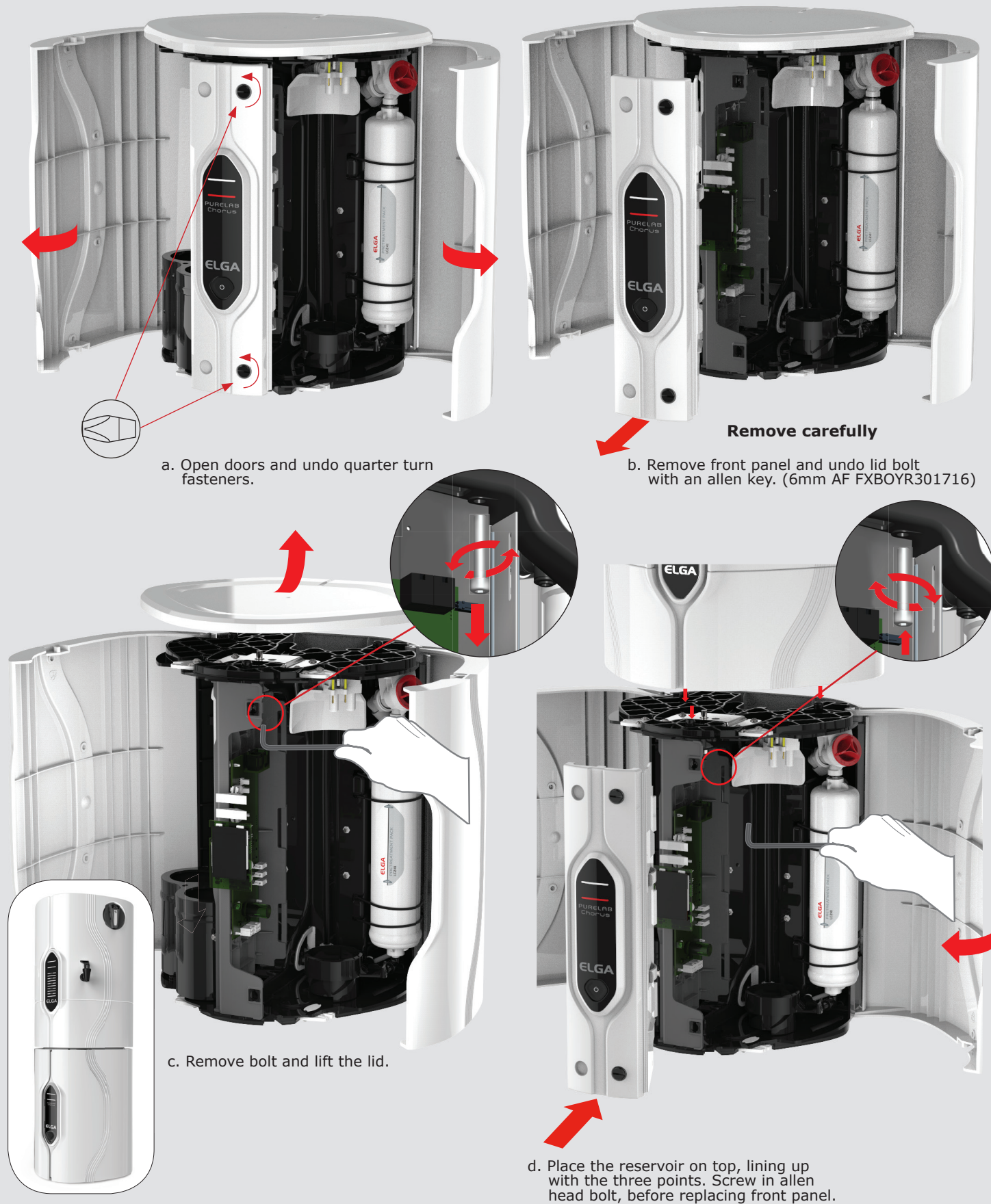
WARNING!

ONLY PLACE ONE CHORUS SYSTEM ON TOP OR BELOW THE RESERVOIR.


SET - UP


1b MODULAR STACKING, RESERVOIR UNDERNEATH (OPTIONAL)

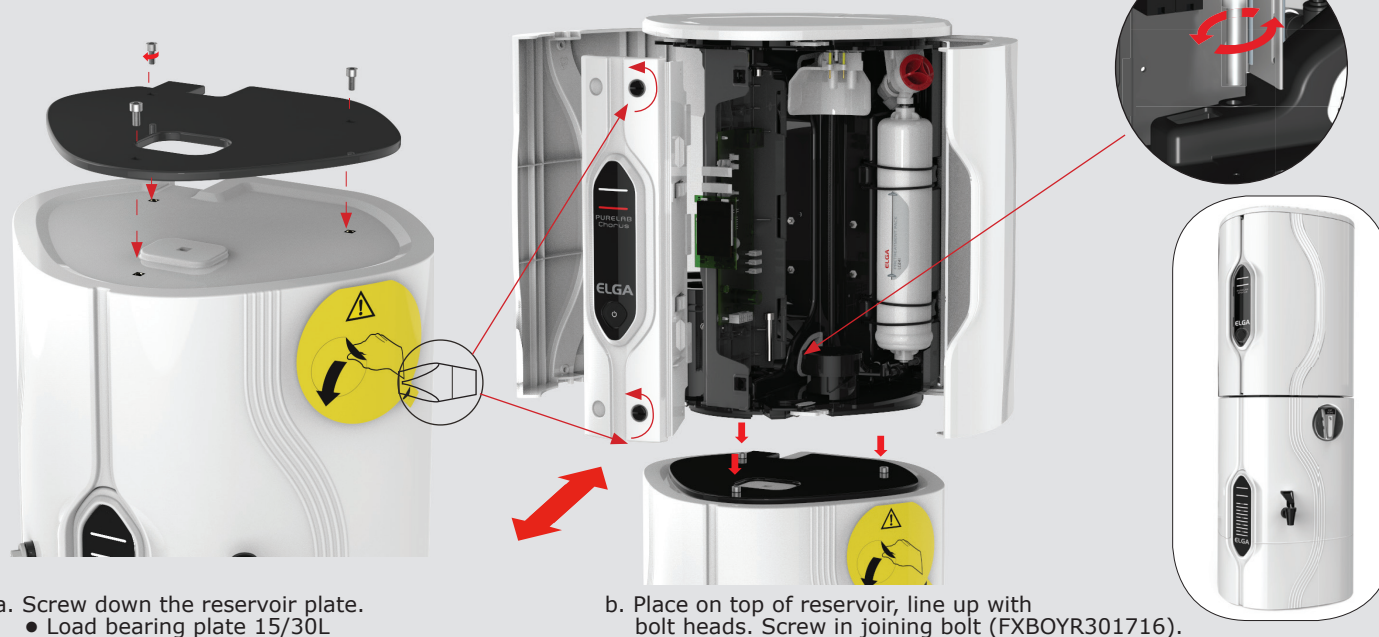2 CONNECT THE WATER (PORT No. 1)

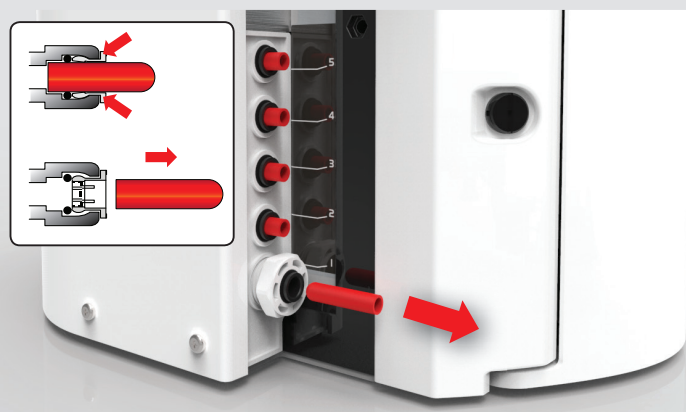
a. Remove the transit plugs from the water inlet connection.

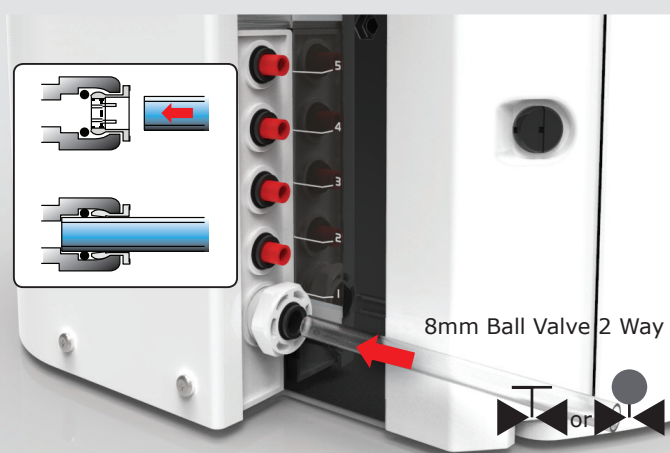
b. Using the tube supplied, firmly push one end securely into the water connector.

3 CONNECT THE DRAIN (PORT No.3)

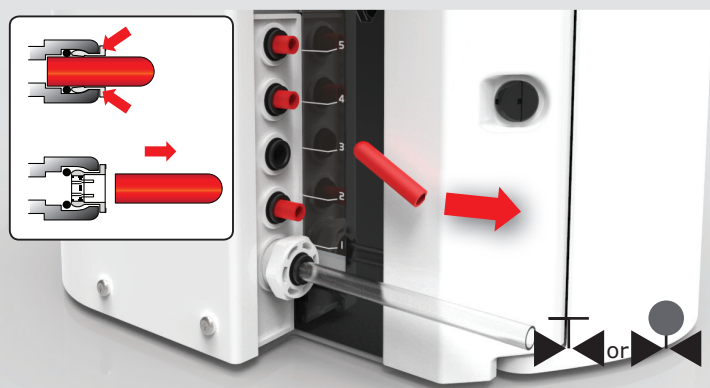a. Remove the transit plugs from the water outlet connection.


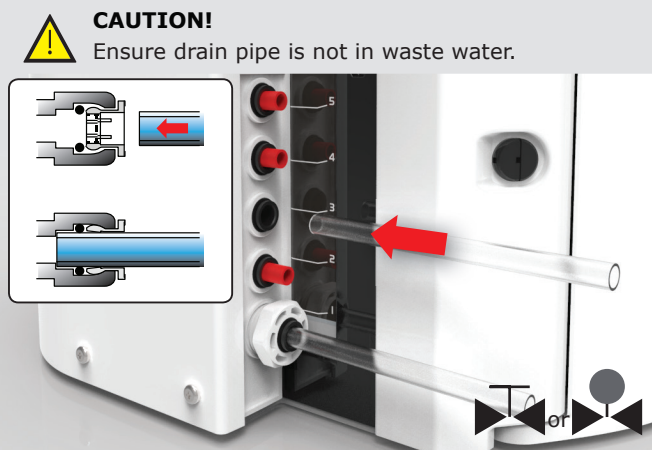
b. Using the tube supplied, firmly push one end securely into the water connector.

Connect the other end of the drain and overflow tubes to a sink or suitable drain capable of handling at least 1.5 l/min. The drain point should have a gravity fall below the level of the unit and any connections direct to drain should have an air-break device fitted.

SET - UP

4 CONNECT THE PURIFIED WATER OUTLET (PORT NO.4)

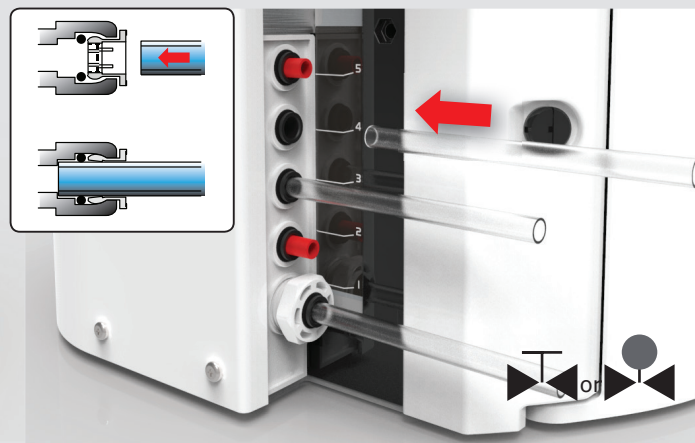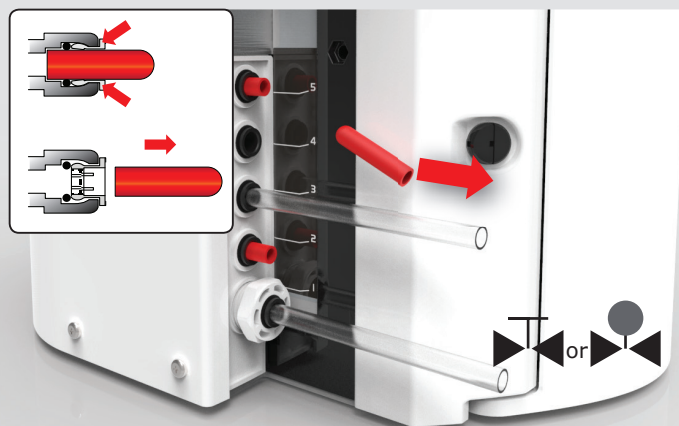a. Remove the transit plugs from the water outlet connection.

b. Using the tube supplied, firmly push one end securely into the water connector.

c. Turn on the water supply.

Minimum Pressure 0.5bar (7.5psi) with boost pump / 4.0bar (58psi) without boost pump.
Maximum Pressure 2.0bar (30psi) with boost pump / 6.0bar (90psi) without boost pump.

N.B This product is not suitable to operate directly with a chorus dispenser.

5 COMMS CONNECTIONS


WARNING!

ONLY USE COMMS CABLE SUPPLIED BY ELGA®.

Connect to:
Reservoir J2 or J1
See 'Reservoir
Installation Guide'.

6 RESERVOIR NOTE

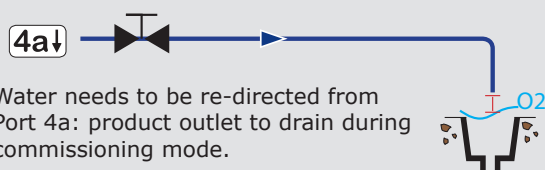
Water needs to be re-directed from Port 4a: product outlet to drain during commissioning mode.

SET - UP

7 REMOVING THE BYPASS PACK LC233

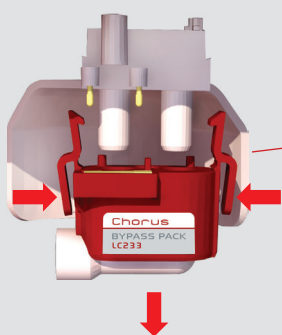Bypass LC233.

a. Open right hand door.

b. Remove Bypass LC233.

c. Remove Purification pack transit plugs

d. Insert Purification pack and turn handle.

e. Make sure it's secure.

8 CONNECT THE POWER (NO.6 PORT)

a. Connect the power cable as shown 1-2-3.


b. Only use the power supply supplied by ELGA®.

The power supply must be earthed.


WARNING!

ONLY USE THE APPLIANCE COUPLER (MAINS LEAD) AND POWER SUPPLY PROVIDED. IF THE EQUIPMENT IS USED IN A MANNER NOT SPECIFIED BY ELGA, THE PROTECTION PROVIDED BY THE EQUIPMENT MAY BE IMPAIRED. POSITION THE POWER SUPPLY SO THAT IT CANNOT COME INTO CONTACT WITH WATER.

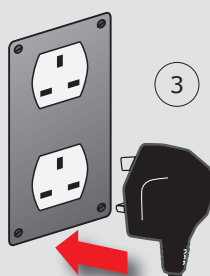3

2

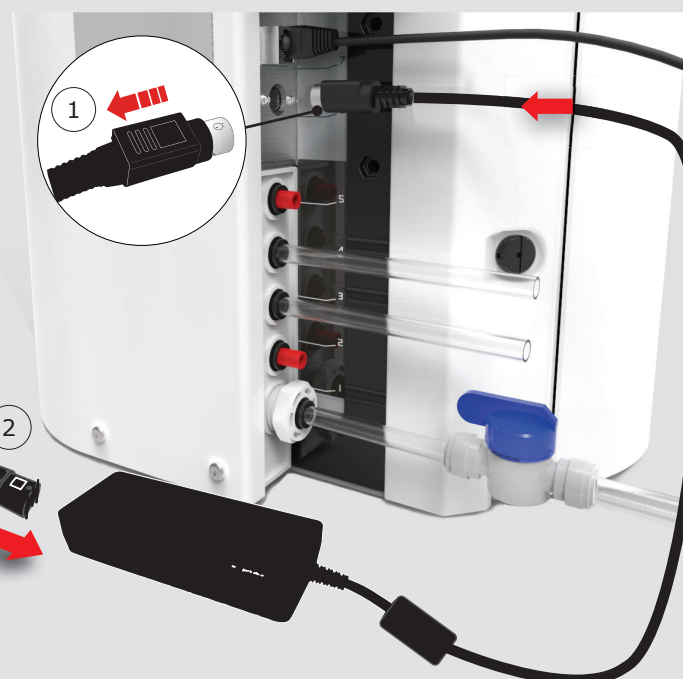1

SET - UP

9 LANGUAGE

a. **Press** for 10 seconds to enter menu.

b. **Press** to scroll, **press and hold** to confirm.


c. **Press once**
White light goes solid.
Allow to operate for approx 60 minutes.


10 FOLLOW THE INSTRUCTIONS

a. After 60 Minutes.

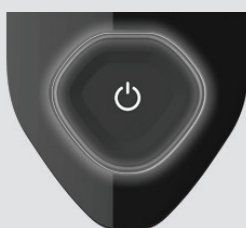b. **Press once**, Process the button (OFF).

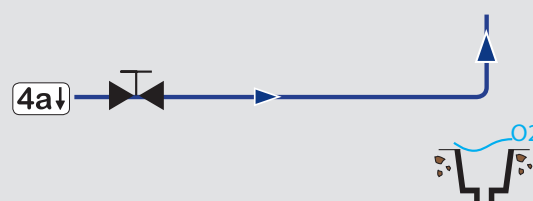c. Reconnect to application or close valve 4a.

d. **Press once**,
White light goes solid.

e. Commissioning is complete and the reservoir will start to fill. The product is now working.

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**TABLE OF CONTENTS**

1.	Introduction	Page 1
1.1	Health & Safety	Page 1
1.2	Product Range	Page 1
1.3	Use of this Manual	Page 1
1.4	Installation	Page 1
1.5	System Configuration	Page 1
1.6	Environment	Page 1
1.7	Commissioning	Page 1
1.8	Customer Support	Page 1
2.	Your Guide to Chorus	Page 2
3.	Maintenance	Page 3
3.1	Cleaning the inlet filter assembly	Page 3
3.2	Replacing the Pre-treatment pack LC241	Page 4
3.3	Replacing the Purification pack LC234	Page 5
3.4	Replacing the RO Reverse osmosis module LC240	Page 6
3.5	Replacing the Composite vent filter LC216	Page 7
4.	Self Help Guidelines	Page 7
5.	Sanitisation Procedure	Page 8 - 10
6.	Health & Safety notes	Page 11
6.1	Environment	Page 11
6.2	Electricity	Page 11
6.3	Residual Biostat	Page 11
6.4	Control Of Substances Hazardous to Health (COSHH)	Page 11
7.	Technical Specifications	Page 12 - 13
8.	Customized Operation	Page 13
8.1	USB Ticket Printing	Page 13
9.	Warranty / Conditions	Page 14

Copyright Note

The information contained in this document is the property of VWS (UK) Ltd, trading as
ELGA® LabWater, and is supplied without liability for errors or omissions.

No part of this document may be reproduced or used except as authorized by contract or other written permission from VWS (UK) Ltd. The copyright and all restrictions on reproducing and use apply to all media in which this information may be placed.

VWS (UK) Ltd. Pursue a policy of continual product improvement and reserve the right to alter without notice the specification, design, price or conditions of supply of any product or service.

© VWS (UK) Ltd. 2019 - All rights reserved.
ELGA® and PURELAB® Are registered trademarks of VWS (UK) Ltd.

Copyright Note

Effersan™ is a registered trademark of Activon Inc.

I. INTRODUCTION

1.1 Health & Safety

Please ensure you read the Health & Safety notes in Section 6.

1.2 Product Range

This Operator manual has been prepared for PURELAB® Chorus 2 RO/DI product models:

- PURELAB® Chorus 2 RO/DI Part No. PC210DIBPM3
- PURELAB® Chorus 2 RO/DI Part No. PC210DIXXM3
- PURELAB® Chorus 2 RO/DI Part No. PC220DIBPM3
- PURELAB® Chorus 2 RO/DI Part No. PC220DIXXM3

1.3 Use of this Manual

This manual guides you through the operation and maintenance of Chorus 2 RO/DI allowing you to obtain a guaranteed supply of purified water to meet your requirements.

1.4 Installation

Quick Start Guide shows you how to install and commission the Chorus 2 RO/DI.

1.5 System Configuration

Maximum system size 4 X Chorus 2 RO/DI + Reservoir.
See Quick Start Guide for typical configurations.

1.6 Environment

The Chorus unit should be installed on a flat, level surface, in a clean, dry environment.
The unit can also be wall mounted using the purpose designed wall mounting kit against a vertical wall capable of supporting the combined weight of the unit and mounting kit INST40989.

**WARNING!**

NOT TO BE USED IN FUME EXTRACTION ENCLOSURES.

1.7 Commissioning



Chorus is supplied in a preset commissioning mode which must be completed before purified water can be dispensed. See Quick Start Guide.

1.8 Customer Support


If you need help with your Chorus, please call your local ELGA representative. For the address of the nearest ELGA LabWater Sales and Service office visit the country list on our website.
www.elgalabwater.com or contact ELGA Lab-Water at: Email: techsupport@elgalabwater.com

2. YOUR GUIDE TO THE CHORUS

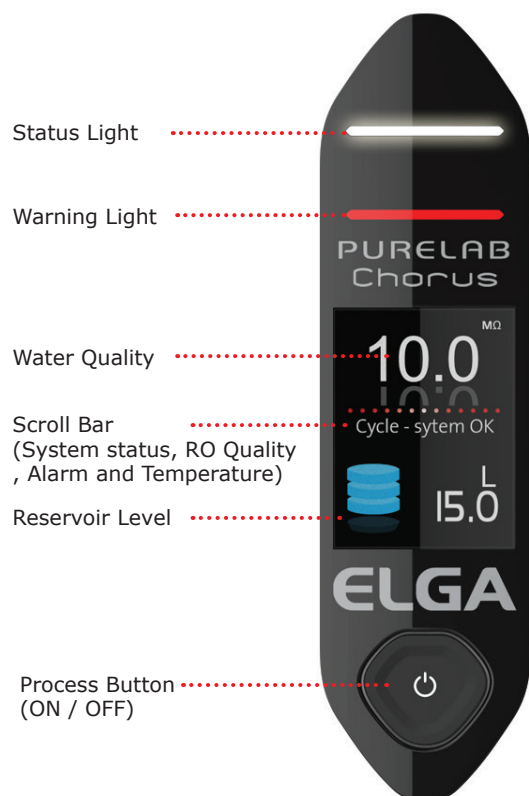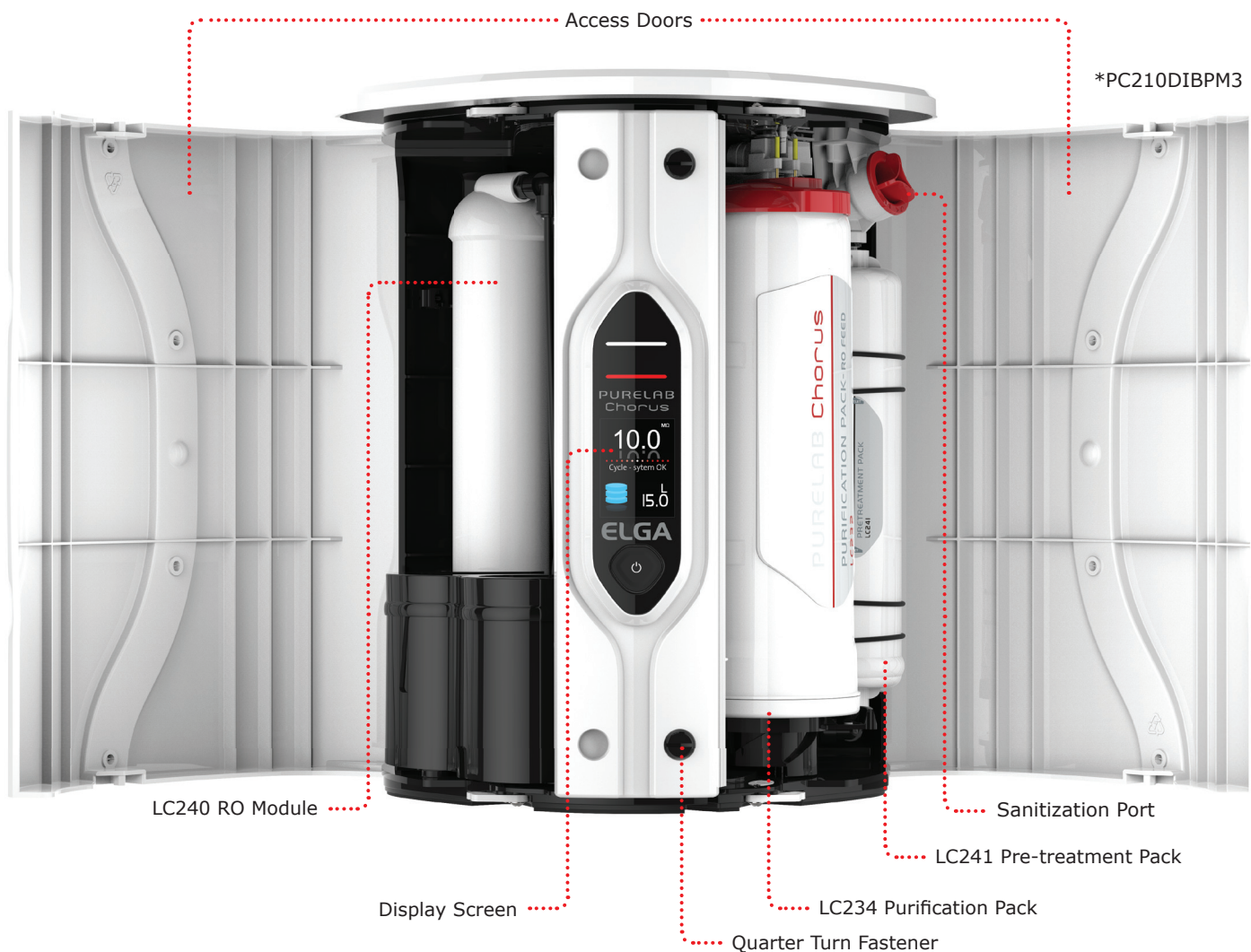1. Water Inlet
2. Drain
3. Product Outlet
4. Return
5. Power In 1
6. COMMS 1
7. Power In 2
8. COMMS 2

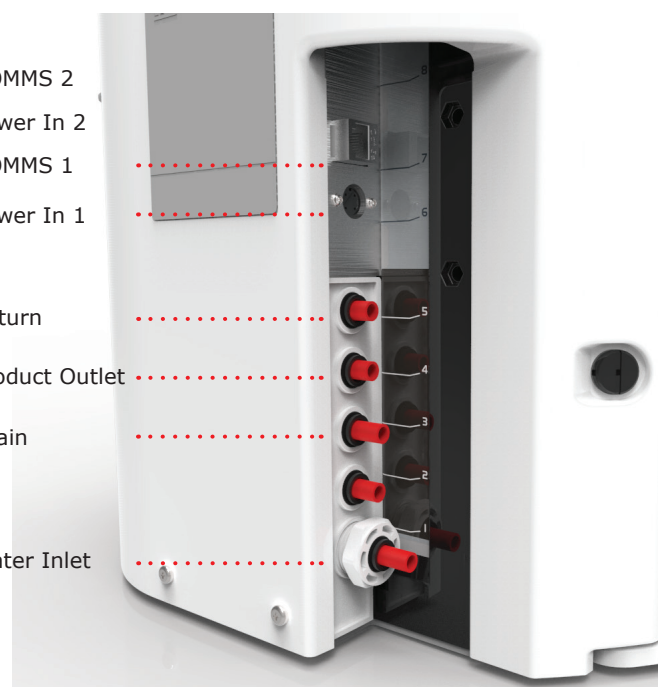3. MAINTENANCE

WARNING!


An approved supplier or distributor should carry out any maintenance work not included in this manual.

Note

Note: Disposal of all end-of-life consumable items should be in accordance with local statutory regulations.

3.1 Cleaning the Inlet Filter Assembly

Follow the instructions below.

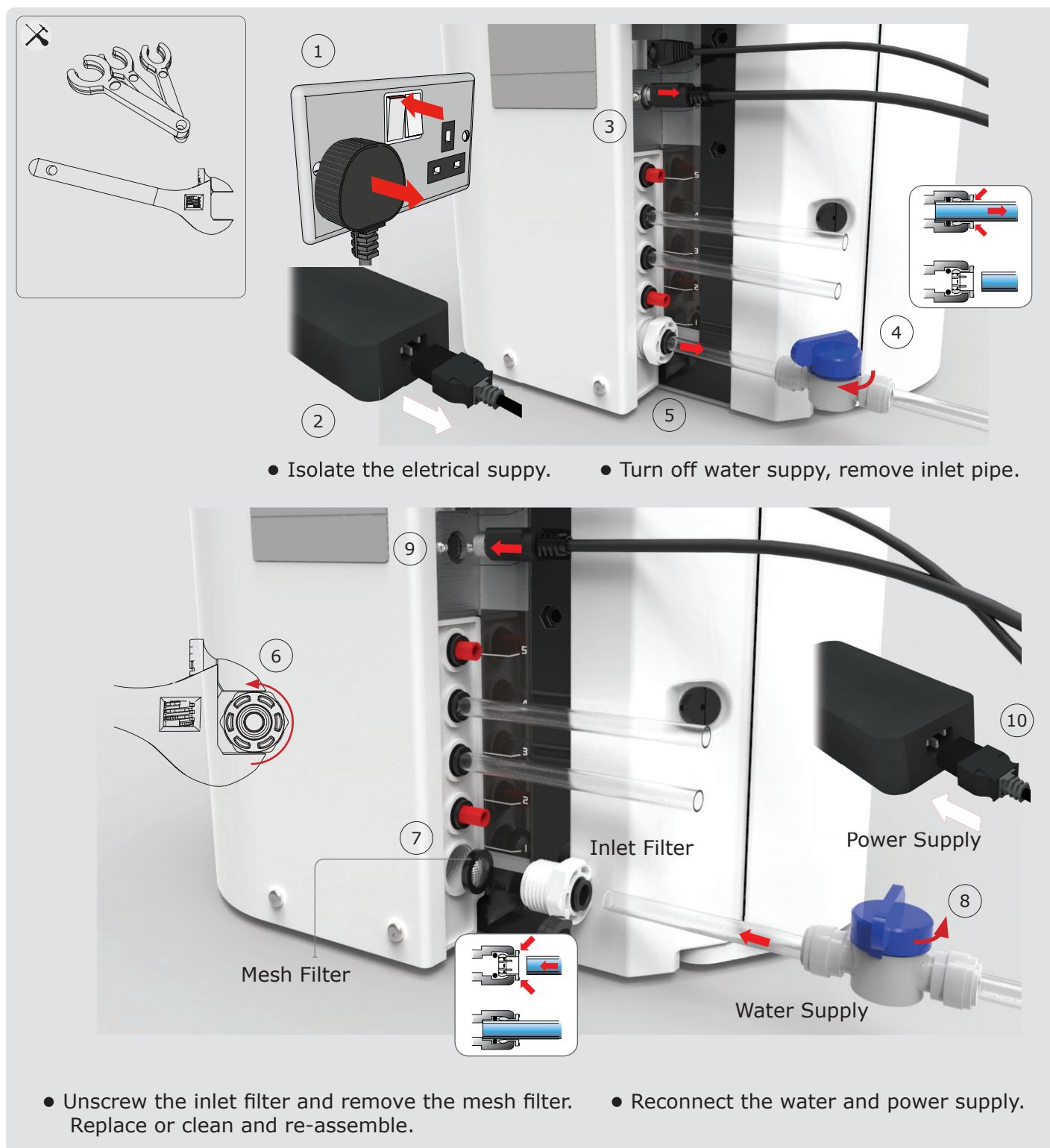1. Isolate the electrical supply.

2. Turn off water supply, remove inlet pipe.

3. Unscrew the inlet filter and remove the mesh filter.

4. Replace or clean and re-assemble.

5. Reconnect the water and power supply.

Mesh Filter

Inlet Filter

Power Supply

Water Supply


3. MAINTENANCE

3.2 Replacing the Pre-treatment Pack LC241

Follow the instructions below.

When 4380 hours (6 mths) is reached the filter change reminder is active. Status light will be **blinking white** showing system needs attention. There is no reminder for replacing the RO Modules (LC240) but it is recommended that they are replaced during this procedure when their service life has expired (2 Yrs).

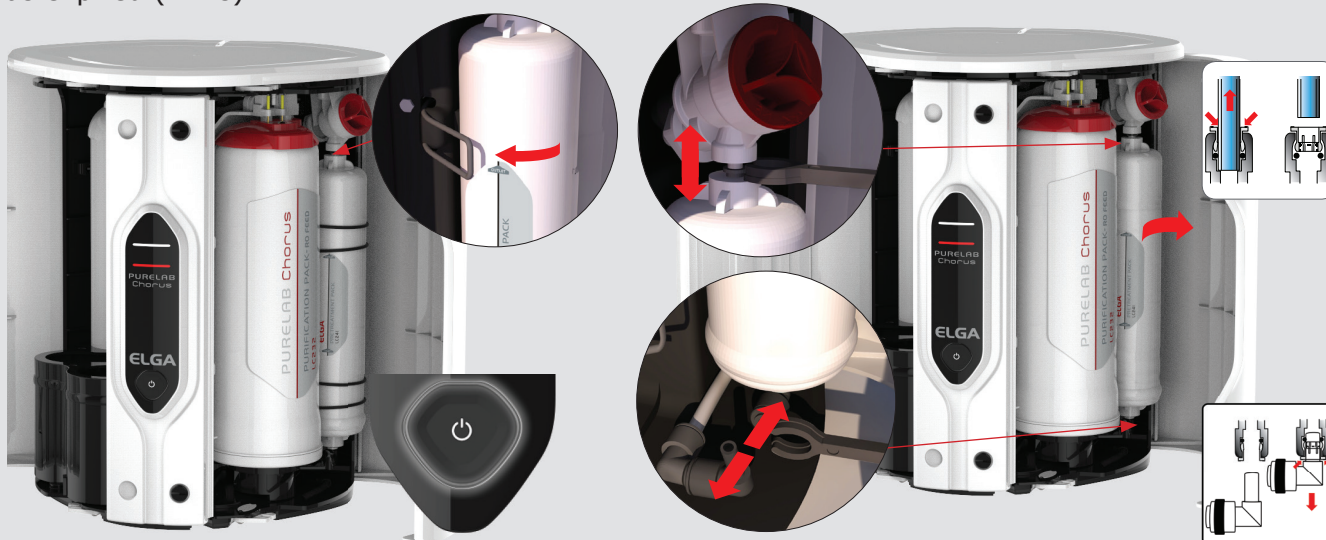1. Press once, Process button (OFF) and unclip Pre-treatment cartridge LC241.


2. Disconnect and remove Pre-treatment cartridge.

3. Remove transit plugs and insert a new Pre-treatment cartridge.

4. Press once, Process button (ON).

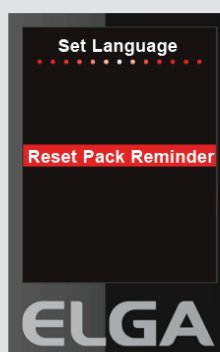5. Press for 10 seconds to enter menu.

6. Press to scroll, press and hold to confirm the reset pack reminder. Normal operation resumes.

3. MAINTENANCE

3.3 Replacing the Purification pack

Follow the instructions below.

When a purification pack is fitted into the system it will record the serial number of the Purification Pack. After 8760 hours (12 months) a reminder is activated to replace the Purification Pack. Status light will be **blinking white** showing system needs attention.


1. Process button (OFF). Rotate the handle to the unlock position, remove Purification Pack.

2. Pull down and tilt towards to remove.


3. Remove transit plugs and insert a new Purification pack.

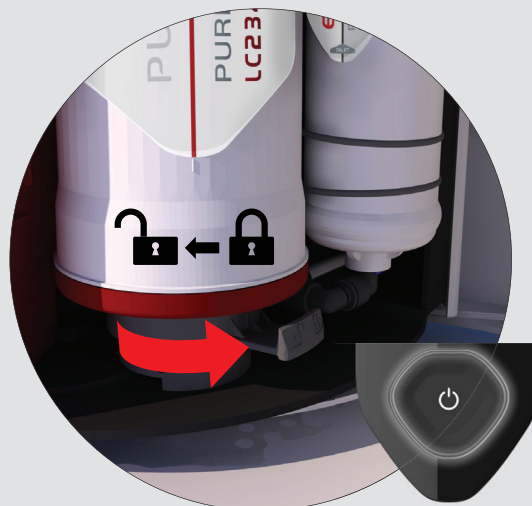4. Process button (ON), normal operation resumes.

5. Press for 10 seconds to enter menu.

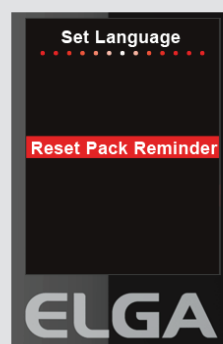6. Press to scroll, press and hold to confirm the reset pack reminder. Normal operation resumes.


3. MAINTENANCE

3.4 Replacing the Reverse Osmosis Module LC240

Follow the instructions below.

Chorus is fitted with Reverse Osmosis (RO) modules. The RO Modules may need replacing if the permeate water purity or flow rate is not adequate and does not meet predicted or previous performance.

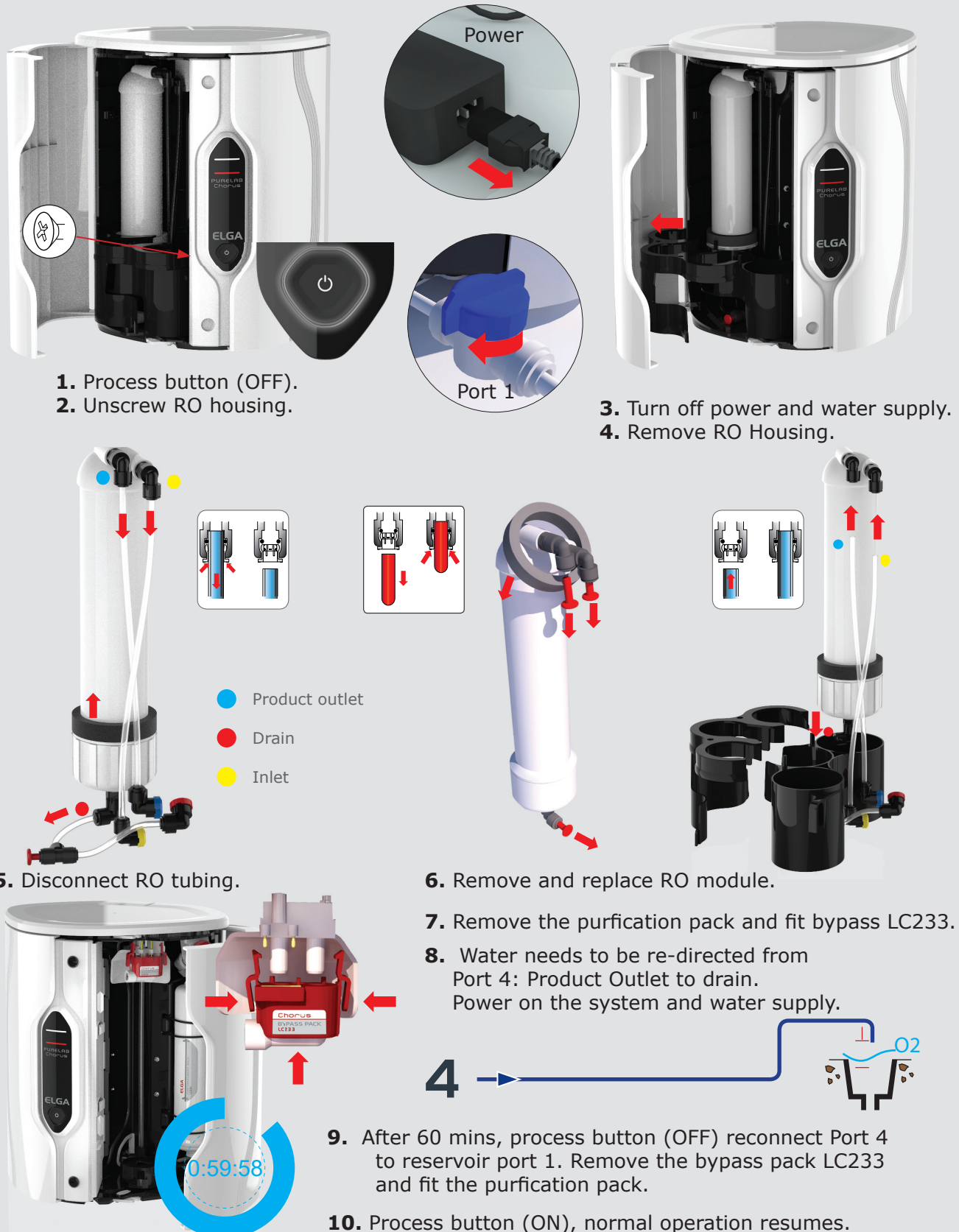WARNING!

ALL NEW RO MODULES ARE FILLED WITH A PRESERVATIVE SOLUTION TO PREVENT BACTERIAL CONTAMINATION DURING STORAGE. THE RO MODULES WILL THEREFORE REQUIRE RINSING PRIOR TO USE.

3. MAINTENANCE

3.5 Replacing the Composite Vent Filter LC216

Follow the instructions below.

4. SELF HELP GUIDELINES

Display Lights

Constant	White	Signifying water purity, unit is operating correctly
Blinking	White	Consumables are nearing the end of useful life
Constant	Red	System needs immediate attention
Blinking	Red / White	USB Software update

Process Button Operation

Press for 10 seconds	Select Languages	Language Menu
Press (At any time)	Start / Stop Modules process	Start and stops modules process
Insert a blank USB and press once. To stop recording data, press once (OFF) before removing the USB.	Load firmware from USB Flash devices	Automatically loads system firmware from USB flash device. Re-start on completion.

Default Alarm Settings

Set Temperature Alarm	Set default values from firmware file. Adjustable set point using application service module. Service module is an engineers tool, only an engineer can adjust the settings. Select to either 30°C, 35°C, 40°C. Default 35°C. https://elgalab.discovertheweb.co.uk
Set Water Purity Alarm Chorus 2	Set default values from firmware file. Adjustable set point using application service module. Resistivity range select from 1.0MΩ.cm, 10.0MΩ.cm or 15.0MΩ.cm. Default Alarm point 15.0MΩ.cm.
Over fill Protection Alarm	Feature activates drain valve if water level is greater than maximum and stops water entering the system to prevent a flood. When alarm is active: V2 energised. If alarm for 5 mins, stop module process. Enter refill cycle on overflow alarm deactivation. Alarm status removed (Reservoir 95%)

5. SANITISATION PROCEDURE

1. Process button (OFF).

2. Water needs to be re-directed from Port 4: Product outlet to drain during sanitisation.

3. Remove purification pack.




4. Change the pre-treatment cartridge LC241.

5. SANITISATION PROCEDURE

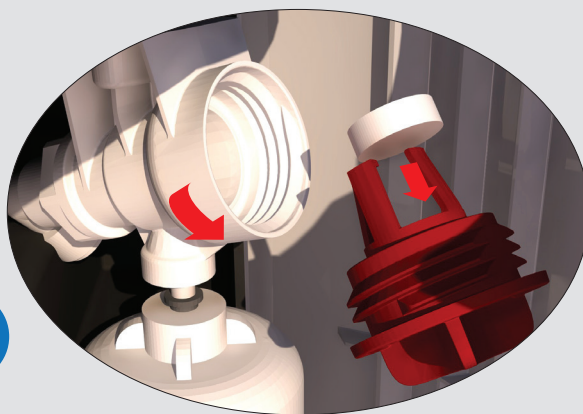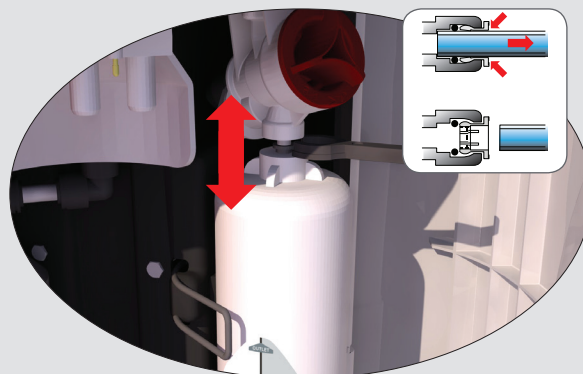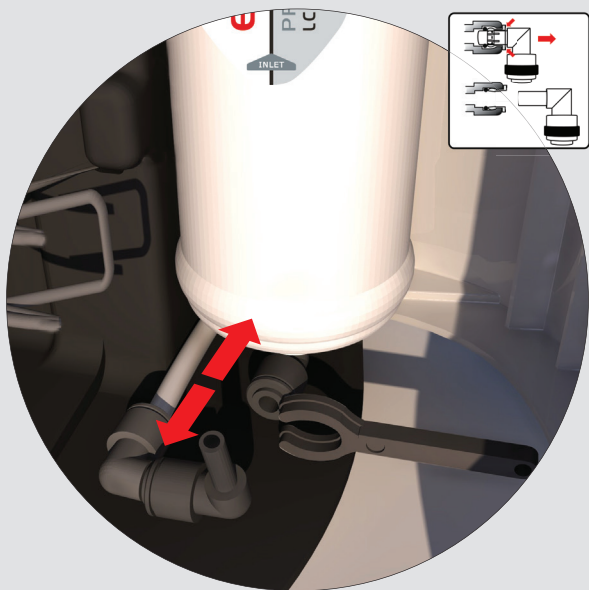
5. Add CT1 or Effersan[™] tablet to the sanitisation port.

Chlorine Tablets:

CT1 Chlorine tablets or EfferSan[™] multi-purpose disinfecting tablets are required to successfully complete this procedure.
(United States EPA Approved EfferSan[™] Tablets)

EfferSan[™] are not supplied and must be sourced from your local ELGA[®] LabWater distributor. Refer to the EfferSan[™] label for additional information and follow all applicable directions for use on the manufacturers label in conjunction with the following instructions.

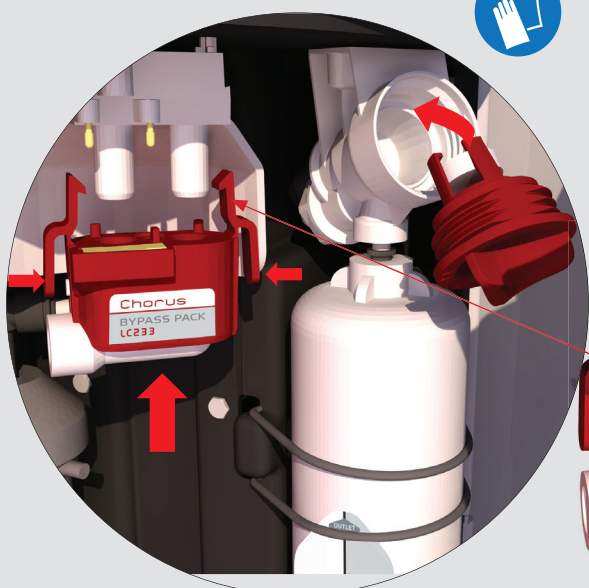6. Fit bypass pack LC233 and Sanitisation cap.

7. Press process button (ON).

5. SANITISATION PROCEDURE

8. Allow to operate for 60 mins.

Note:

*Please be aware that during this procedure the pre-treatment quality alarm will be active.

Port 4

4

Reservoir Port 1

9. After 60 mins, press the process button (OFF) and reconnect port 4 to reservoir port 1.

10. After 60mins remove bypass pack LC233 and insert a purification pack.

11. Process button (ON), normal operation will now resume.

6. HEALTH & SAFETY NOTES

**WARNING!**

WARNINGS ARE GIVEN WHERE FAILURE TO OBSERVE THE INSTRUCTIONS COULD RESULT IN INJURY OR FATALITY.

**CAUTION!**

Cautions are given where failure to observe the instructions could result in damage to the equipment, associated equipment and processes.

6.1 Environment

**CAUTION!**

Clean dry indoor. Temp 5-40°C. Humidity max 80% non-condensing.
Failure to follow the environmental specification could result in damage to the system.
The Chorus is not designed for use in fume cupboards where chemicals could damage the system.

**WARNING!**

ALL NEW RO MODULES ARE FILLED WITH A PRESERVATIVE SOLUTION TO PREVENT BACTERIAL CONTAMINATION DURING STORAGE.
THE RO MODULES WILL THEREFORE REQUIRE RINSING PRIOR TO USE.

6.2 Electricity

The appliance coupler (mains lead) or power supply connected to the rear of the unit can be removed to isolate the power supply. If access to this is restricted then it is recommended that access to the supply socket is easily available to disconnect the electrical supply.

**WARNING!**

ONLY USE THE APPLIANCE COUPLER (MAINS LEAD) AND POWER SUPPLY PROVIDED.
THE USE OF THESE WILL ENSURE ADEQUATE EARTH PROTECTION IS PROVIDED.

IF THE EQUIPMENT IS USED IN A MANNER NOT SPECIFIED BY ELGA®, THE PROTECTION PROVIDED BY THE EQUIPMENT MAY BE IMPAIRED.

POSITION THE POWER SUPPLY SO THAT IT CANNOT COME INTO CONTACT WITH WATER

**WARNING!**

ALWAYS ENSURE THAT THE ELECTRICAL POWER SUPPLY IS ISOLATED BEFORE WORKING INSIDE CHORUS.

6.3 Residual Preservative

**WARNING!**

DURING THE COMMISSIONING CYCLE TRACE RESIDUAL PRESERVATIVE IS FLUSHED FROM THE SYSTEM.
A MATERIAL DATA SHEET IS AVAILABLE UPON REQUEST.

6.4 Control of Substance Hazardous to Health (COSHH)

Material safety data sheets are available upon request.

7. TECHNICAL SPECIFICATIONS

Purified Water Specification					
Model	PURELAB® Chorus 2 PC210DIXXM3	PURELAB® Chorus 2 PC210DIBPM3	PURELAB® Chorus 2 PC220DIXXM3	PURELAB® Chorus 2 PC220DIBPM3	
Make-up rate @15°C (*4bar feed)	10l/h	10l/h	20l/h	20l/h	
Delivery flow rate - maximum	240l/24h	240l/24h	480l/24h	480l/24h	
Inorganic Rejection	1 to >10MΩ. cm	1 to >10MΩ. cm	1 to >10MΩ. cm	1 to >10MΩ. cm	
Organics (TOC) - Typical	<30ppb	<30ppb	<30ppb	<30ppb	
Organics (MW >200 Dalton)	>99% Rejection	>99% Rejection	>99% Rejection	>99% Rejection	
Bacteria - Typical	<5 CFU/ml	<5 CFU/ml	<5 CFU/ml	<5 CFU/ml	
Particles	>99% Rejection	>99% Rejection	>99% Rejection	>99% Rejection	
Feed Water Specification					
Source Quality	Potable mains supply	Potable mains supply	Potable mains supply	Potable mains supply	
Fouling Index-Maximum	10	10	10	10	
TDS-Maximum	1000ppm	1000ppm	1000ppm	1000ppm	
Conductivity	<2000µS/cm	<2000µS/cm	<2000µS/cm	<2000µS/cm	
Free/total Chlorine-Maximum	0.5ppm	0.5ppm	0.5ppm	0.5ppm	
Chloramine - Maximum	0.5ppm	0.5ppm	0.5ppm	0.5ppm	
Heavy Metals-Maximum (Fe/Mn)	0.05ppm	0.05ppm	0.05ppm	0.05ppm	
Silica-Maximum	30ppm	30ppm	30ppm	30ppm	
Temperature	1-35°C	1-35°C	1-35°C	1-35°C	
Flowrate (max requirement @ 15°C)	100l/h	100l/h	100l/h	100l/h	
Drain requirements	80l/h	80l/h	80l/h	80l/h	
Maximum feed pressure	6.0 bar (90 psi)	2.0bar (30psi)	6.0 bar (90 psi)	2.0bar (30psi)	
Minimum feed pressure	4.0 bar (58 psi)	0.5bar (7.5psi)	4.0 bar (58 psi)	0.5bar (7.5psi)	

Connections	
Inlet	8mm (5/16) OD tube
Drain	8mm (5/16) OD tube
Outlet	8mm (5/16) OD tube

7. TECHNICAL SPECIFICATIONS

Electrical Requirements					
Model	PURELAB® Chorus 2 PC210DIXXM3	PURELAB® Chorus 2 PC210DIBPM3	PURELAB® Chorus 2 PC220DIXXM3	PURELAB® Chorus 2 PC220DIBPM3	
Mains Input	100-240 Vac 50-60 Hz	100-240 Vac 50-60 Hz	100-240 Vac 50-60 Hz	100-240 Vac 50-60 Hz	
System control voltage	24V dc	24V dc	24V dc	24V dc	
Power consumption (peak demand)	60VA	60VA	60VA	60VA	
Reservoir Level Connection	RJ45	RJ45	RJ45	RJ45	
Noise level	<45dBA	<45dBA	<45dBA	<45dBA	
Dimensions & Weight					
Model	PURELAB® Chorus 2 PC210DIXXM3	PURELAB® Chorus 2 PC210DIBPM3	PURELAB® Chorus 2 PC220DIXXM3	PURELAB® Chorus 2 PC220DIBPM3	
Height	435mm	435mm	435mm	435mm	
Width	375mm	375mm	375mm	375mm	
Depth	340mm	340mm	340mm	340mm	
Supply Weight	17kg(37lb)	20kg(44lb)	17kg(37lb)	20kg(44lb)	
Operating Weight	22kg(49lb)	22kg(49lb)	22kg(49lb)	22kg(49lb)	
Installation	Bench/Wall/Floor	Bench/Wall/Floor	Bench/Wall/Floor	Bench/Wall/Floor	

* Dependent on feed water as part of our policy of continual improvement we reserve the right to alter the specifications given in this document.

Consumables			
Part No.	Description	Typical Service Life *	Max. Shelf Life
LC241	Pre-treatment Pack	6 Months (4380 hours)	2 years
LC240	RO Module	2 years	2 years
LC234	Purification Pack	12 Months (8760 hours)	2 years
LC216	Composite Vent Filter	6 Months	2 years
CT1	Tablets	1 every 6 months	2 years
Effersan™	Tablets	1 every 6 months	>8 Months

*Service life is an estimate only, and will depend on the application and feed water quality. Please take care to ensure your order is the correct consumables items.

8. CUSTOMIZED OPERATION

It's possible to customize your Chorus by downloading additional features and changing your current settings. To find out more, please visit: <https://elgalab.discovertheweb.co.uk>

8.1 USB Data logging To download technical data please follow instructions.

CAUTION!
USE A BLANK USB.

1. Insert a blank USB and press the process button.
2. To stop recording data, press process button (OFF) before removing the USB.

9. WARRANTY / CONDITIONS

ELGA® LabWater is a trading name of VWS (UK) Ltd. General Limited Warranty

VWS (UK) Ltd. warrants the products manufactured by it against defects in materials and workmanship when used in accordance with applicable instructions for a period of one year from the date of shipment for the products. VWS (UK) Ltd. MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED. THERE IS NO WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.

The warranty provided herein and the data, specifications and descriptions of the VWS (UK) Ltd. products appearing in VWS (UK) Ltd. published catalogues and product literature may not be altered except by express written agreement signed by an officer of VWS (UK) Ltd. Representations, oral or written, which are inconsistent with this warranty or such publications are not authorized and, if given, should not be relied upon.

In the event of a breach of the foregoing warranty, VWS (UK) Ltd. sole obligation shall be to repair or replace, at its option, any product or part thereof that proves to be defective in materials or workmanship within the warranty period, provided the customer notifies VWS (UK) Ltd. promptly of any such defect. The exclusive remedy provided herein shall not be deemed to have failed of its essential purpose so long as VWS (UK) Ltd. is willing and able to repair or replace any nonconforming VWS (UK) Ltd. product or part. VWS (UK) shall not be liable for consequential, incidental, special or any other indirect damages resulting from economic loss or property damage sustained by any customer from the use of its products.

VWS (UK) Ltd. Warranty

VWS (UK) Ltd. warrants the water systems manufactured by it, BUT EXCLUDING MEMBRANES AND PURIFICATION PACKS, against defects in materials and workmanship when used in accordance with the applicable instructions and within the operating conditions specified for the systems for a period of one year from the earlier of:

the date of installation, or the 120th day following the date of shipment.

VWS (UK) LTD. MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED. THERE IS NO WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. The warranty provided herein and the data, specifications and descriptions of the VWS (UK) Ltd. systems appearing in VWS (UK) Ltd. published catalogues and product literature may not be altered except by express written agreement signed by an officer VWS (UK) Ltd. Representations, oral or written, which are inconsistent with this warranty or such publications are not authorised and, if given, should not be relied upon.

In the event of a breach of the foregoing warranty, VWS (UK) Ltd. sole obligation shall be to repair or replace, at its option, any product or part thereof that proves to be defective in materials or workmanship within the warranty period, provided the customer notifies VWS (UK) Ltd. promptly of any such defect. The cost of labour for the first ninety (90) days of the above warranty period is included in the warranty; thereafter, labour cost shall be at the customer's expense. The exclusive remedy provided herein shall not be deemed to have failed of its essential purpose so long as VWS (UK) Ltd. is willing and able to repair or replace any nonconforming VWS (UK) Ltd. system or component part. VWS (UK) Ltd. shall not be liable for consequential, incidental, special or any other indirect damages resulting from economic loss or property damage sustained by any customer from the use of its process systems.

Products or components manufactured by companies other than VWS (UK) Ltd. or its affiliates

("Non - VWS (UK) Ltd. products") are covered by the warranty, if any, extended by the Product manufacturer. VWS (UK) Ltd. hereby assigns to the purchaser any such warranty; however VWS (UK) LTD. EXPRESSLY DISCLAIMS ANY WARRANTY WHETHER EXPRESSED OR IMPLIED, THAT THE NON - VWS (UK) Ltd. PRODUCTS ARE MERCHANTABILITY OR FIT FOR A PARTICULAR PURPOSE.

NOTICE

VWS (UK) Ltd. is constantly striving to improve its products and services. Consequently, the information in this document is subject to change without notice and should not be construed as a commitment by VWS (UK) Ltd.

Also VWS (UK) Ltd. assumes no responsibility for any errors that may appear in this document. This manual is believed to be complete and accurate at the time of publication. In no event shall VWS (UK) Ltd. be liable for incidental or consequential damages in connection with or arising from the use of this manual.

VWS (UK) Ltd. warrants its products against defects in materials and workmanship as described in the Warranty statement.

The Labwater Specialists

ELGA is the global laboratory water brand name of Veolia. The information contained in this document is the property of VWS (UK) Ltd, trading as ELGA LabWater, and is supplied without liability for errors or omissions. © VWS (UK) Ltd. 2019-All right reserved. ELGA®, PURELAB®, ADEPT® Labpure® and PureSure® are registered trademarks of VWS (UK) Ltd.

Tel: +44 203 567 7300 Email: info@elgalabwater.com Web: www.elgalabwater.com

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© 2019 ELGA® Labwater / VWS (UK) Ltd. All rights reserved.