

ISS.	MN/PRN/DN	ACTION	APPR'D	DATE
1	DN/E/144	JP		

CONNECTION & DIMENSION DETAILS.

REF	FUNCTION	DETAILS	LA586 25 LITRE	LA587 40 LITRE	LA590 75 LITRE
A	INLET	8mm OD TUBE CONNECTOR	✓	✓	✓
B	OUTLET	8mm OD TUBE CONNECTOR	✓	✓	✓
C	OVERFLOW	8mm OD TUBE CONNECTOR	✓	✓	✓
D	DISPENSE TAP	—	✓	✓	✓
E	LEVEL CONTROL	DIN PLUG / LEAD - 6 WAY	✓	✓	✓
F	RESERVOIR ACCESS	80mm DIAMETER	✓	✓	✓
OPTIONAL WALL MOUNTING KIT - ELGA PART No.			LA591	LA591	LA592
OVERALL HEIGHT - MM			500mm	600mm	755mm
DIMENSION 'A'			415mm	525mm	680mm
DIMENSION 'B'			360mm	400mm	510mm
DIMENSION 'C'			270mm	300mm	355mm
OVERALL HEIGHT - INCHES			19.68"	23.60"	29.70"
DIMENSION 'A'			16.33"	20.60"	26.77"
DIMENSION 'B'			14.17"	15.70"	20.07"
DIMENSION 'C'			10.60"	11.80"	13.97"

LEVEL CONNECTION DETAILS


DETAILS	HIGH LEVEL CONNECTION		INTERMEDIATE LEVEL CONNECTION		LOW LEVEL CONNECTION	
	N/O	COM	N/C	COM	N/C	COM
PURELAB PRIMA 60/120	BLACK	RED	WHITE	YELLOW	-	-
PURELAB OPTION-S 30/60	BLACK	RED	WHITE	YELLOW	-	-
PURELAB OPTION-R 30/60	BLACK	RED	WHITE	YELLOW	BLUE	GREEN
PURELAB OPTION-E 25/50/80	BLACK	RED	WHITE	YELLOW	-	-
MEDICA 30/60/100	BLACK	RED	WHITE	YELLOW	-	-
MEDICA-R 30/60/100	BLACK	RED	WHITE	YELLOW	-	-
PURELAB ULTRA RANGE	-	-	-	-	-	-

NOTE:
THE BLUE AND GREEN WIRES ARE NOT USED IN MANY CASES AND SHOULD BE INSULATED FROM EACH OTHER.

UNIT WEIGHTS

UNIT	WEIGHT (EMPTY)	WEIGHT (FULL)
LA586 (25 LITRE)	5.0KG (11.0LB)	30.0KG (66.0LB)
LA587 (40 LITRE)	8.0KG (17.6LB)	48.0KG (105.6LB)
LA590 (75 LITRE)	12.0KG (26.4LB)	87.0KG (191.0LB)

IF IN DOUBT-ASK

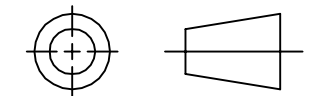

3rd ANGLE PROJECTION DO NOT SCALE

MATERIAL: N/A

FINISH: N/A

GENERAL NOTE: N/A

TOLERANCES UNLESS STATED:
ONE PLACE DEC. ± N/A ANGULAR ± 0.5°
TWO PLACE DEC. ± ALL HOLE TOLS. TO BS.4500 H11

DIMS. IN mm	SCALE	DRAWN	JP	AUG 11	DATE
SIZE: A3	NTS	CHECKED	SB	AUG 11	
		APPROVED	LU	AUG 11	

TITLE
**SALES DIAGRAM
OUTLINE & CONNECT'N
LA586/LA587/LA590
3 SWITCH RESERVOIRS**

PART NUMBER: N/A


ELGA

www.elgalabwater.com


info@elgalabwater.com

ELGA is the laboratory water brand of Veolia Water Solutions & Technologies.


DRG. NO.	SD_3SR_01	ISSUE	1
SHEET No.1 OF 1			

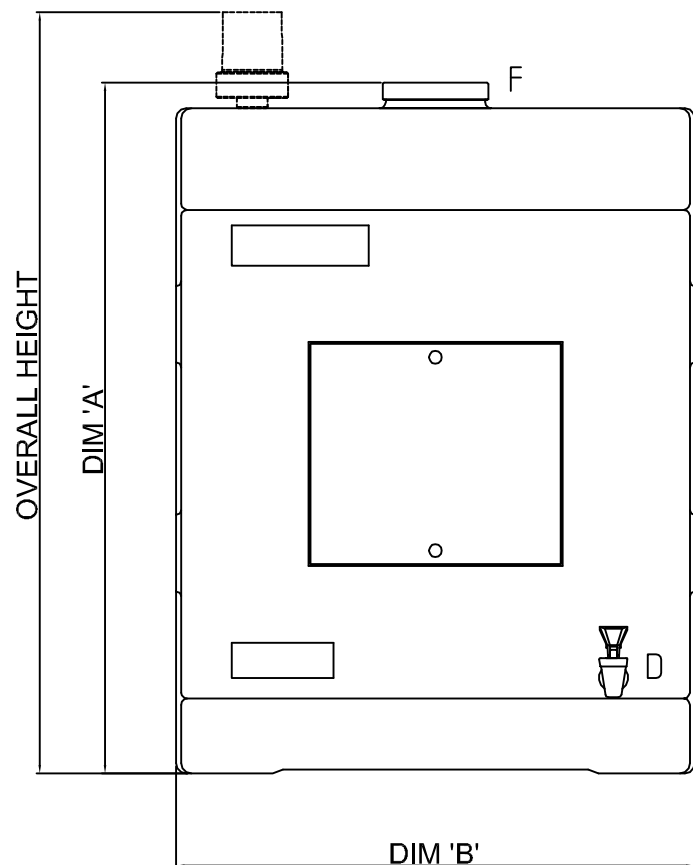

FRONT VIEW

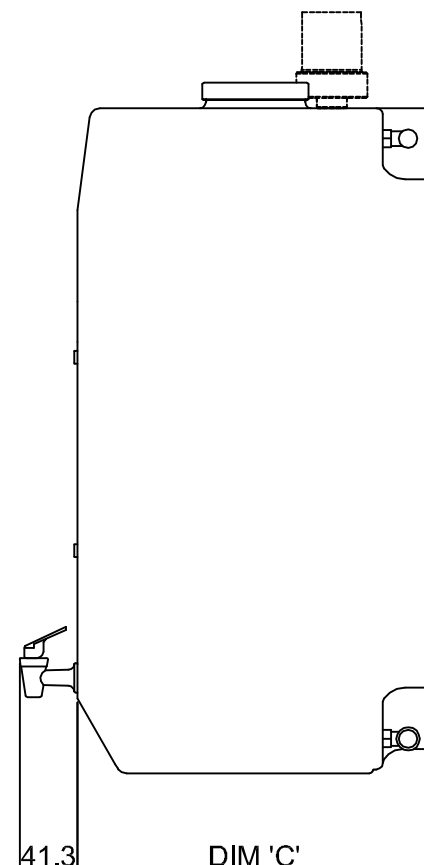

SIDE ELEVATION

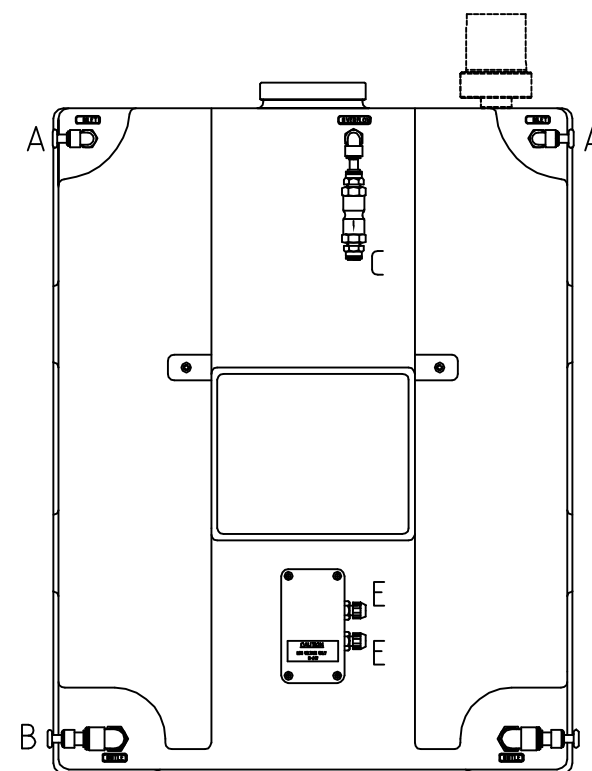

REAR VIEW

THIS DRAWING IS THE PROPERTY OF ELGA Labwater. THE DETAILS SHOWN ARE STRICTLY CONFIDENTIAL AND MUST NOT BE DISCLOSED TO ANY THIRD PARTY WITHOUT THE CONSENT OF VWS (UK) LTD.